

Operating instructions

Model-No.: JK 162 / 52-4 TT HP

1019102-00H / us-en / 01.2022



Contact information



Customer Service JK Products & Services

411 West Washington Avenue
Suite B
Jonesboro, Arkansas 72401
800-445-0624
870-935-9046
<http://jkproducts.us/service>

Manufacturer JK-Products GmbH

Rottbitzer Str. 69
53604 Bad Honnef
GERMANY
Tel.: +49 (0) 22 24 / 818-140
Fax: +49 (0) 22 24 / 818-166

WARNING!



Read these operating instructions.

Failure to observe these instructions:

- may result in serious injury and death,
- can result in damage to the device and the environment.
- For the safety of yourself and others, read and observe all instructions contained herein before operating the device.
- Please make the instructions and additional information available to anyone who operates or services the device.

Copyright

The copyright remains with JK-Holding GmbH.

The contents may not be duplicated (in whole or part), distributed or utilized for competition purposes without authorization or disclosed to other parties.

The right to make technical changes in regard to the figures and information contained in this instruction manual is reserved!

1 Table of Contents



2	What certain words and signs mean (Glossary).....	7
2.1	Definitions.....	7
2.2	Warning signs.....	8
2.3	Symbols on the device.....	10
3	What this device is.....	12
3.1	Guidelines.....	12
3.2	Export.....	12
3.3	Scope of delivery.....	12
3.3.1	Optional	12
3.4	Device description.....	13
3.5	Options	14
3.6	Function description.....	15
3.7	Transport	15
3.8	Installation location	16
3.8.1	Ventilation.....	16
3.8.2	Setup location.....	17
3.8.3	Commissioning.....	18
3.9	Signs and stickers on the device.....	19
3.10	Technical data	30
3.10.1	Power rating, electrical specifications and noise level.....	30
3.10.2	Dimensions.....	32
3.10.3	Lamp configuration	33
3.10.4	Accent lighting	34
3.10.5	Spare parts.....	34
4	Why you use it.....	35
4.1	Intended use.....	35
4.2	Skin condition	35
5	Who cannot use it	36
5.1	Prohibited usage	36
5.1.1	Tanning system.....	36
5.1.2	Vibra Shape (vibration plate).....	36
6	What you must do to avoid serious harm.....	37
6.1	General	37
6.2	Tanning system.....	38
6.3	Sensor.....	39



7	What you must do to avoid other harm.....	40
7.1	General.....	40
7.2	Tanning system.....	40
7.3	Sensor	40
7.4	Vibra Shape (vibration plate).....	40
7.5	Foreseeable misuse	41
7.6	Obligations of the studio operator.....	41
7.7	Staff training	41
8	Risks of using this device	42
9	Benefits of using this device	43
10	How to decide to use this device	44
11	What to do before the treatment	45
12	How to use this device.....	46
12.1	How to set this device up	46
12.2	How to make sure the device works	46
12.3	How to use it to treat your condition	46
12.4	How to choose correct treatment regimen	46
12.4.1	Skin types.....	46
12.4.2	Tanning times.....	47
12.5	How to operate the device.....	48
12.5.1	Control Center.....	48
12.5.2	User settings overview	49
12.6	Functions.....	51
12.6.1	Navigation	53
12.6.2	Description	54
12.7	Selection of the application	56
12.8	Operating display of the sensor	57
12.9	Sun Select setting with the outer interface	57
12.10	Measurement (automatic mode)	57
12.11	Connections	58
12.11.1	Couple with Bluetooth® devices	58
12.11.2	MP3 player and smartphone	60
12.12	How to use	61
12.12.1	How to start	61
12.12.2	Application tips	61
12.12.3	How to stop	61
12.13	What to do after the treatment	62
12.13.1	Cleaning and disinfection	62
12.14	Safety instructions for cleaning/disinfection and maintenance	63
12.15	How to clean/disinfect.....	66
12.15.1	Cleaning and disinfection of surfaces	67
12.15.2	Cleaning/disinfection plan	68



12.16	How to maintain.....	69
12.16.1	Technical modifications	69
12.16.2	Maintenance plan.....	69
12.16.3	Cleaning and maintenance work preparation.....	72
12.16.4	Instructions for replacing the lamps	77
12.16.5	Cleaning and disinfection of the sensor.....	82
12.16.6	Cleaning or replacing lamps in the base	82
12.16.7	Cleaning or replacing the lamps in the side section	83
12.16.8	Cleaning or replacing lamps and filter panels in the canopy	84
12.16.9	Cleaning or replacing the lamps and filter panels in the shoulder tanner	86
12.16.10	Cleaning or replacing the interior lights	88
12.16.11	Cleaning or replacing the canopy accent lighting	88
12.16.12	Cleaning or replacing the front panel accent lighting.....	88
12.16.13	Cleaning the filters in the base	88
12.16.14	Cleaning the canopy filter mats.....	89
12.16.15	Cleaning the filter mats in shoulder tanner.....	89
12.16.16	Air conditioning: Empty the condensation tank and clean the filter	90
12.16.17	AQUA / AROMA SYSTEM maintenance.....	93
12.16.18	Cleaning the air intake openings.....	95
12.17	How to store the device	96
12.17.1	Decommissioning	96
12.17.2	Storage	96
12.17.3	Disposal	96
13	What to do when things go wrong	97
13.1	Troubleshooting.....	99
14	Where you can find out more.....	100
15	How we studied this device.....	101
16	Appendix	102
16.1	Model-No.: JK 162 / 52-4 TT HP: EMC information.....	102
16.1.1	Table 201: Guidelines and manufacturer's declaration – Electromagnetic emissions.....	102
16.1.2	Table 202: Guidelines and manufacturer's declaration – Electromagnetic interference.....	103
16.1.3	Table 204: Guidelines and manufacturer's declaration – Electromagnetic interference.....	104
16.1.4	Table 206: Recommended safe distances between portable and mobile HF communications equipment and Model-No.: JK 162 / 52-4 TT HP	105



17	JK-Products, Inc. – Limited Warranty	106
17.1	Limited Acrylic Warranty.....	107
17.2	Warranty Claims	108
17.2.1	Policy & Procedures	108
17.2.2	Labor Warranty Claim/Credit Procedures	109
17.2.3	Quality inspection and acceptance form	109
17.3	Limited Lifetime Warranty.....	110
18	Index	111

2 What certain words and signs mean (Glossary)



2.1 Definitions

Contraindications

Issues where the risk is always greater than the benefit and there are data supporting the risk and benefit relationship.

Studio operator

People who provide tanning devices to users on a commercial basis. The studio operator is trained to operate and maintain the device and is responsible for the orderly operation of the device as well as compliance with the maintenance intervals.

User

A person who uses the tanning device at a commercial location.

Staff

People who are responsible for the operation, cleaning/disinfection as well as general maintenance work and who instruct the users in the operation of the devices.

JK Factory Certified Technician

Level 1

Authorized to perform general maintenance such as replacing lamps, aqua mist, and fragrance containers as well as overall cleaning of the unit.

Level 2

Authorized to perform maintenance and repairs such as replacement of ballast, framed acrylic, electronic circuit boards, and other electrical components.

Level 3

Authorized and trained to utilize JK Products proprietary diagnostic software to adjust operating parameters and perform detailed diagnostics.

Level 4

This individual is uniquely qualified to assess the condition of a unit and make a determination if the device meets the requirements for continued use ten (10) years from the date of manufacture.



NOTE:

Each technician level is authorized to perform the tasks of technician levels smaller in number to theirs. E.g. a level 3 technician is authorized to perform the tasks assigned to a level 1 and 2 technician.

2.2 Warning signs



Danger notices, nonobservance will lead to death or serious injury:

DANGER!



Type and source of danger, e.g. 'DANGER – High voltage!'

Description of the consequences, e.g. 'You can be killed or injured by an electric shock!'

Countermeasures:

- Disconnect device from power supply before opening.
- Do not reach into the inside with the device switched on.

Warning notices, nonobservance can lead to death or serious injury.

WARNING!



Type and source of danger, e.g. 'WARNING – hot surface!'

Description of the consequences, e.g. 'You can burn your hand (skin, face)!'

Countermeasures:

- Do not touch hot parts and nearby surfaces.
- Wait until hot parts and surfaces have cooled down.

Warning notices, ultraviolet radiation, nonobservance can lead to serious injury:

WARNING!



Type and source of danger, e.g. 'WARNING – ultraviolet radiation!'

Description of the consequences, e.g. 'You can burn your skin!'

Countermeasures, e.g.:

- Avoid overexposure.

WARNING!



Crushing hazard, e.g. when closing the bed base!

Fingers can be jammed.

- Do not hold fingers between bed base and front panel.



Caution notices, nonobservance can lead to less than serious injury or damage of property.

⚠ CAUTION!



Unauthorized lamps!

The word 'Caution' together with the warning sign points out that less than serious injury can be expected here.

Countermeasures:

- Avoid doing this.

NOTICE



ESD-sensitive components!

The word 'NOTICE' without warning sign points out that danger for devices, materials and the environment can be expected here.

Countermeasures:

- Avoid doing this.

Important information:



NOTE:

This symbol does not identify safety notices, but states information for a better understanding of the processes.

2.3 Symbols on the device



The following symbols are displayed on the device or on the packaging:



Manufacturer, company address



Type BF applied part¹

Qty:

Quantity



Article Number/Reference Number



Serial Number

+ UDI

Unique Device Identification



Healthcare Barcode



Read the operating instructions!




Therefore you must never fill the AQUA SYSTEM canister with tap water or other liquids.



Do not add any liquid that remains in the old canister to the new AQUA SYSTEM canister!



Note the expiration date on the AQUA SYSTEM canister.

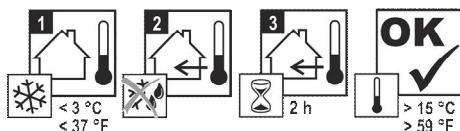

Where there are great differences in temperature between the transport route and the installation site the device may not be put into operation immediately after being set up.

Wait at least 2 hours before connecting the device to the power supply.

Maximum transport weight



Do not stack!



¹ Devices with physical contact between device and user belong to category BF (body float).



This side up.



Keep dry.



Temperature limitations.

The temperature limits are noted next to the upper and lower horizontal line.



Recyclable material.



Polystyrene



Cardboard



Mixed paper



Steel



Wood



Permitted air pressure during transport and storage.

3 What this device is



3.1 Guidelines

This product is in conformity with the performance standards for sunlamp products under 21 CFR 1040.20, FDA and listed with ETL to UL 482 and CAN/CSA 22.2 No. 224.

3.2 Export

This device is exclusively intended for use in the USA and Central America, and may not be exported to other countries and operated there! The manufacturer assumes no liability whatever in the event of non-compliance. High liability risks may occur for the exporter and/or the studio operator in the event of noncompliance.


3.3 Scope of delivery

- Tanning device
- Lamps
- Shoulder tanner with audio system (with Voice Guide)
- Enclosed tanning tunnel
- Service kit including an Allen-key and suction cup for filter panels
- Air conditioning system
- Protective goggles
- Mastercard NFC (studio operator)
- Technical documentation (folder with instructions, brochures for error codes and default settings and additional documentation)

3.3.1 Optional

- Central exhaust air socket


3.4 Device description

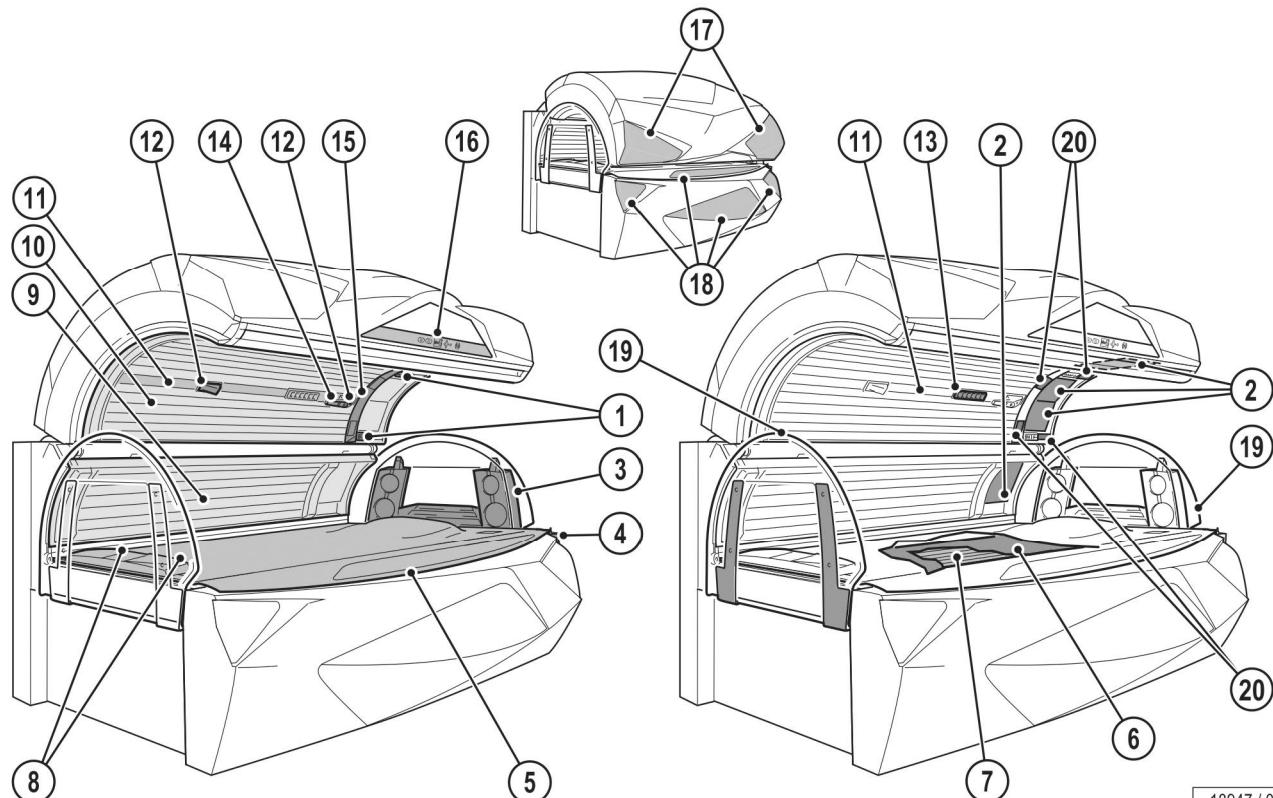

18947 / 0

1. Adjustable air nozzles – Facial cooling
2. Facial tanner (high pressure lamps + low pressure lamps)
3. Shoulder tanner with audio system
4. Headphones and AUX connection
5. Acrylic panel base
6. Partition
7. Low pressure lamps, base
8. Body cooling air nozzles, foot end
9. Low pressure lamps, side section
10. Low pressure lamps, canopy
11. Interior lighting (LED)
12. Nozzles AQUA / AROMA SYSTEM
13. Body cooling air nozzles
14. Operating handle
15. Display unit
16. Control Center
17. Accent lighting, canopy (LED)
18. Accent lighting, front panel (LED)
19. Enclosed tanning tunnel (optional)
20. Hybrid Performance LED

3.5 Options

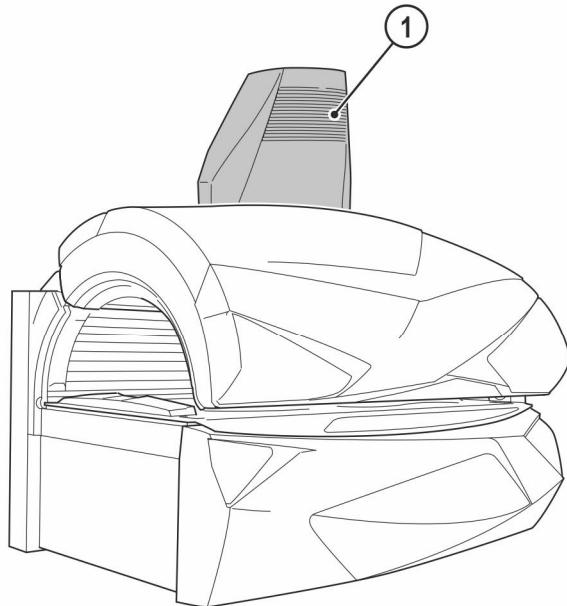

14732 / 0

1. Central exhaust air socket

Exhaust air (Option)

As an optional comfort feature, your device can be equipped with an exhaust duct which allows the warm air from the device to be ducted outside of the room or building. Doing so can be an effective way in some climates to reduce the energy costs associated with Air conditioning while at the same time increasing not only user comfort, but the environmental comfort of the entire studio as well.

Things to consider before purchasing the exhaust option are:

1. Climate; extremely hot climates are typically not a good choice for exhausting air outside due to the need to cool down the hot air coming in to replace the air being exhausted, however mild to cooler climates can benefit from exhausting the warm air from the device thereby allowing cooler air to enter from outside the building.
2. Device location; in some installations, walls or other obstruction may prevent the use of an exhaust duct.



NOTE:

- Always consult a licensed heating, ventilation, and air conditioning (HVAC) Technician when considering the option to exhaust the warm air from your device. They can help you decide if an exhaust is best for your specific location.
- The exhaust is a comfort option and does not increase or decrease the safety and effectiveness of the device. The exhaust duct is not required to keep the device in a certain temperature range in order to operate safe and effectively.



3.6 Function description

Tanning system

The primary technical components of a tanning device are an artificial source of UV radiation, a variety of filters and reflectors as well as a mechanical structure with a defined active surface.

Different tanning results can be achieved in varying strength tanning devices. This is due to the different strength UV lamps and the different UV-A and UV-B proportions of the UV radiation. The UV-A proportion primarily generates a superficial tan, which appears rapidly and is intensive but also fades more rapidly, the UV-B radiation is primarily responsible for more long-term tanning results.

Traditional UV light with the addition of (beauty) red light results in a natural looking tan. Optimal tanning results are achieved due to the synergistic combination of light spectrums.


NFC CONNECT¹

Near field communication (NFC) is a communication standard between smartphones and similar devices that allows contactless data exchange. Communication is also possible between an NFC device and an unpowered NFC chip.

The NFC point of this tanning device for example may be used to store basic settings on a chip card. These settings can be read for the next session and may be changed if necessary.

3.7 Transport

Tanning system

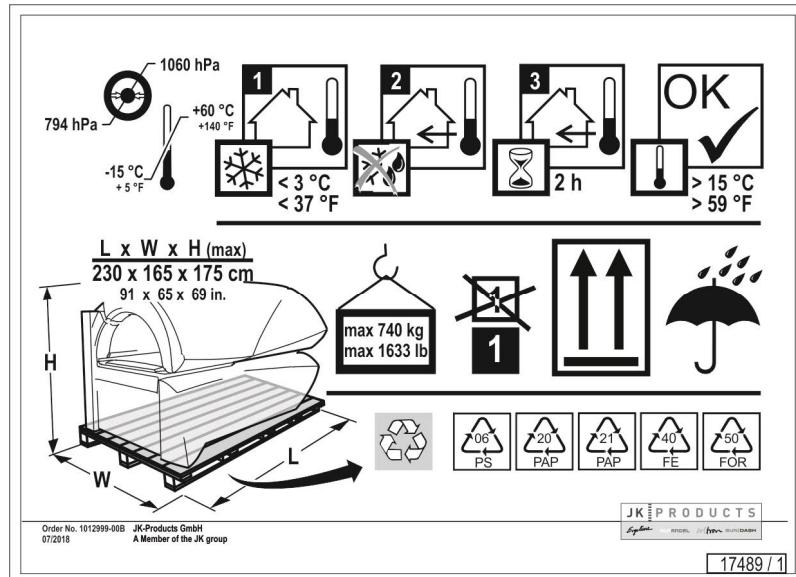

¹ Communication system: NFC (Near Field Communication); ISO/IEC 14443A/MIFARE; max. operating range: approx. 50 mm, frequency range: 13.56 MHz



3.8 Installation location

This sunlamp product is suitable for tanning salons, fitness studios or similar installation locations.

You may find detailed information on installation and setting up in the Installation Instructions (see folder Technical Documentation).


3.8.1 Ventilation

NOTICE

Danger of failure!

Danger of overheating by obstruction of the air routing.

- Observe the ventilation information to avoid serious system malfunctions.

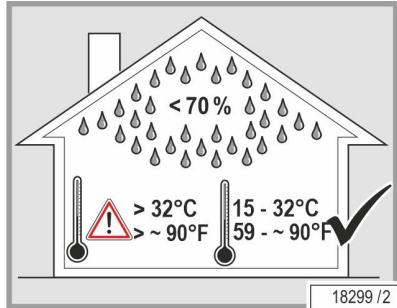

You can obtain data regarding the air supply and ventilation from the customer service department (see page 2).



3.8.2 Setup location

⚠ WARNING!



Sudden movements and malfunctions!

Crushing hazard.

- The device is not suitable for operation in mobile facilities (ships, buses, trains).

⚠ CAUTION!



Danger of overheating by obstruction of the air routing!

Any liability is excluded in case of damage or injury proven to be the result of the use of non-genuine spare parts.

Use of non-compliant lamps voids warranty and risks injury!

- Do not change, obstruct or block the air inflow and outflow area to the device or make any unauthorized changes to the device. The manufacturer shall bear no liability for damage or injuries resulting from such actions.
- Observe the ventilation information to avoid serious system malfunctions.
- Before putting a device into operation, the transport pallet must be removed from under the device. Danger of overheating by obstruction of the air routing.
- Maintain the minimum distances from the walls! Dimensions are given on page 32.
- Only use OEM spare parts of the same type! The use of other parts will void the manufacturer's warranty!

NOTICE

Salty air! Contaminated air!


Device damage due to corrosion of the housing and electrical components.

- Do not install the device near swimming pools.
- Do not operate the device in locations with aggressive environmental conditions (e.g. air containing chlorine).

Damage to the device possible!

Where there are great differences in temperature between the transport route and the installation site the device may not be put into operation immediately after being set up.

- Wait at least 2 hours before connecting the device to the power supply.

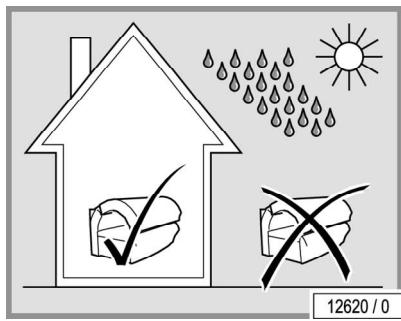

NOTICE

Damage to the device possible!
The device is not to be used outdoors.

Strong Radio Frequency (RF) emitters!

Please do not install or operate the device in close proximity to strong Radio Frequency (RF) emitters such as metal detectors or electromagnetic anti-theft systems. Please inform the installation personal of any existing RF emitting systems prior to installation.



NOTE:

Sources of RF emissions are not always obvious, therefore special care should be taken to ensure all sources have been properly identified before installation of the device.

Environmental conditions

- Optimum ambient temperature:
between 59 °F (15 °C) and ~ 90°F (~ 32°C)
- Storage temperature: between 5 °F (-15 °C) and 140 °F (+60 °C)
- Air pressure: between 794 hPa and 1060 hPa
- Relative humidity: between 30 % and 70 %

Please contact the customer service department before putting the device into operation - see page 2.

3.8.3 Commissioning

DANGER!



Damaged mains power cable!

Danger to persons through electric shock or burns.

- Only the manufacturer or specialist companies are permitted to replace the mains power cable.

Initial commissioning will be undertaken by the manufacturer's own service personnel or by a (Level 3) JK Factory Certified Technician. The device will be handed over ready for use.

As the studio operator you are responsible for abiding by the electrical regulations applicable at the installation site.



The unit incorporates special preventive measures with regard to the electromagnetic compatibility and must be installed and commissioned according to the EMC instructions obtained chapter 16 Appendix.

If the device has not been operated or serviced for a period of one year or longer, it must first be inspected by a (Level 3) JK Factory Certified Technician before being put back into operation. Failure to do so may result in personal injury or property damage.




3.9 Signs and stickers on the device

Signs which warn of danger areas or give important information on components are attached to the device. Ensure that the instruction signs are always clearly recognizable and easily legible. Any missing instruction signs or stickers must be replaced.

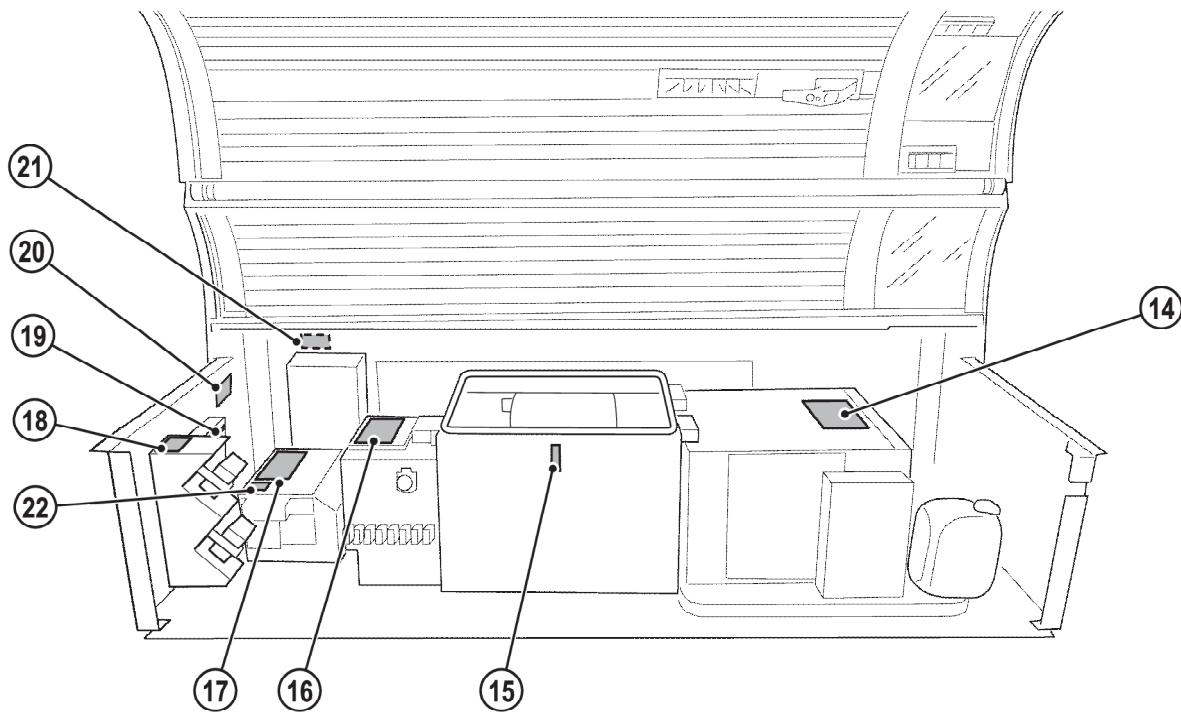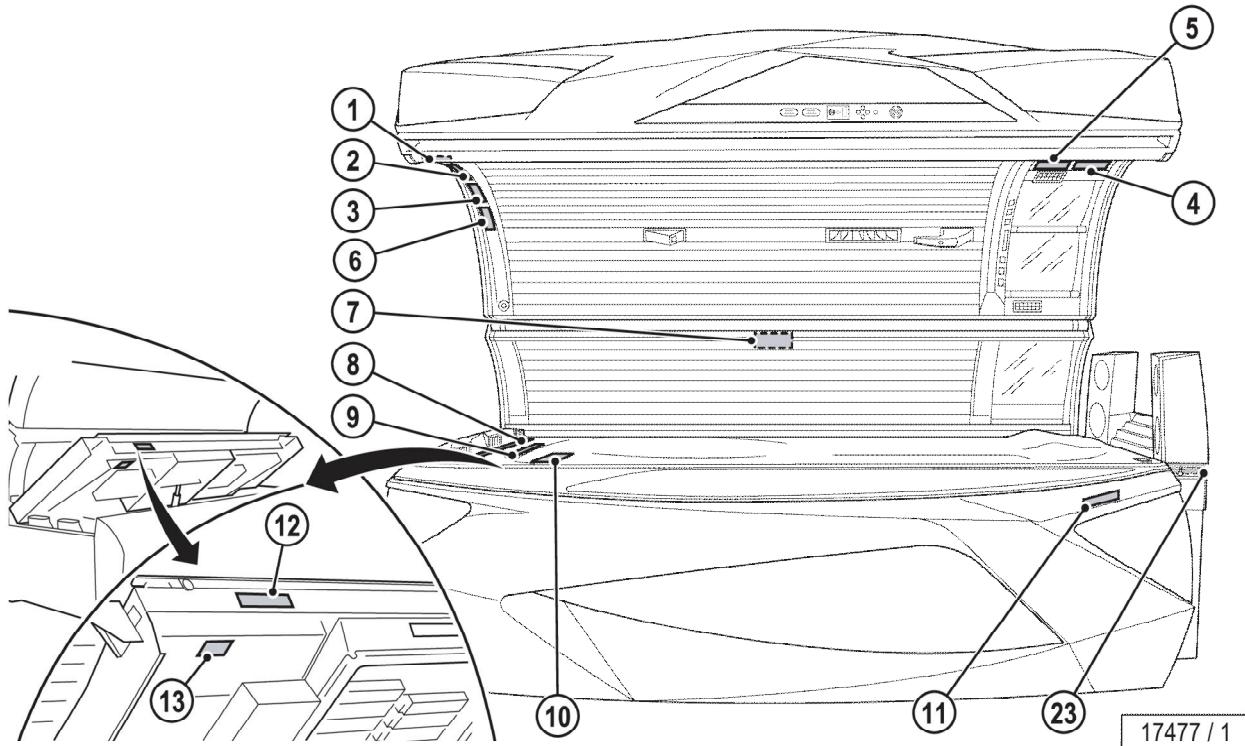




1: Name plate (example)

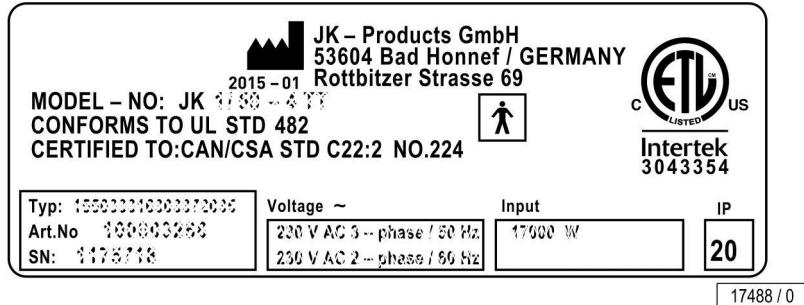

The name plate is attached to the interior side of the canopy. It contains important information for identifying the device (e. g. the serial number = Serial-No.) and the conformity mark.

2: UDI plate (example)

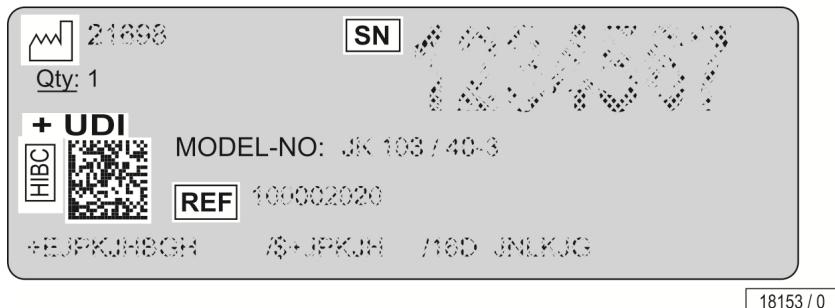

3: Standards




4: Sticker 'Follow Operating Instructions' (1014232-..)





5: Warning sticker UV radiation (801140-..)

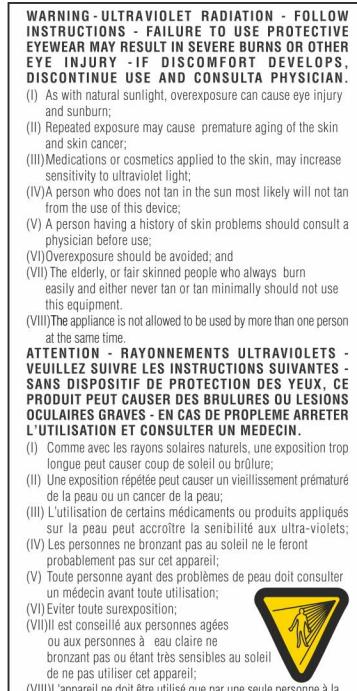


12076 / 2

6: Sticker (1014367-..)



17197 / 2

7: Sticker (843319-..)

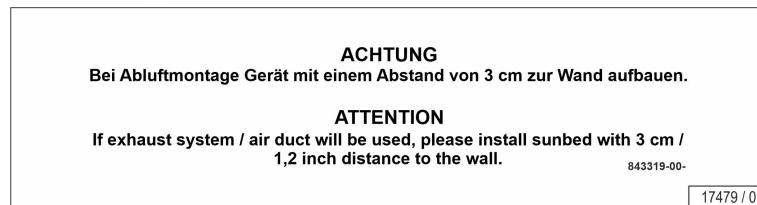

17479 / 0

The sticker is mounted on the rear panel.

8: Exposure label 1019217-..

DANGER ULTRAVIOLET RADIATION		Recommended Exposure Schedule					
Skin Type	Week 1 1st - 3rd Treatments	Week 2 4th - 6th Treatments	Week 3 7th - 10th Treatments	Week 4 11th - 15th Treatments	Weekly Subsequent Treatments		
I - Sensitive	Not Advised	Not Advised	Not Advised	Not Advised	Not Advised	Not Advised	
II - Fair	2 Min.	4 Min.	6 Min.	8 Min.	10 Min.		
III - Average	2 Min.	4 Min.	6 Min.	8 Min.	10 Min.		
IV - Brown	2 Min.	6 Min.	8 Min.	10 Min.	10 Min.		
V - Dark Brown	2 Min.	6 Min.	8 Min.	10 Min.	10 Min.		

WEAR PROTECTIVE EYEWEAR! Failure to do so may result in severe burns or long term injury to the eyes.	This unit utilizes UV lamps. Replace only with Ergoline ULTRA 520W (high-pressure / facial tanner) and Ergoline ULTRA 250 W (high-pressure / shoulder tanner) and Genesis Type U HYBRID Performance 200W or other recognized equivalents.
Certain medications or cosmetics may increase your sensitivity to the ultraviolet radiation. Consult a physician before using sunlamps if you are using medication, have a history of skin problems or believe yourself to be especially sensitive to sunlight. If you do not tan in the sun, you are unlikely to tan from the use of this product.	To obtain the recommended exposure distance, lie in the middle of the acrylic glass panel and lower the canopy to the closed position. Canopy will stop at the correct distance from lounge acrylic surface. Warning: The use of any other position may result in overexposure!

Tanning can begin on a regular basis. A tan normally appears after several exposures and is maximized after four (4) weeks of exposure in accordance with the recommended exposure schedule for your skin type. In order to maintain your tan 10 minute sessions are recommended from the fifth week onward.	Use protective eyewear whenever unit is energized. Read and follow instruction booklet before using this product. To avoid electrical shock during maintenance disconnect power from unit. This product is in conformity with performance standards for sunlamp products under 21CFR1040.20, FDA.
Do not tan more than once in a 24-hour period.	

Part No: 1019217XX

18534 / 2

The sticker is mounted on the lamp cover foot end.

9: Lamp cover base sticker (801139-..)

CAUTION - RISK OF PERSONAL INJURY UV LIGHT SOURCE - KEEP PROTECTIVE BARRIER IN PLACE	ATTENTION - RISQUE DE BLESSURE SOURCE DE RAYONS UV - NE PAS ENLEVER L'Ecran
Keep protective barrier in place or replace removable parts after servicing. Ne pas enlever l'écran protecteur ou remettre en place les pièces amovibles après l'entretien.	


801139-00a

17192 / 0

The sticker is mounted on the lamp cover foot end.



10: Sticker partition (1002317-..)

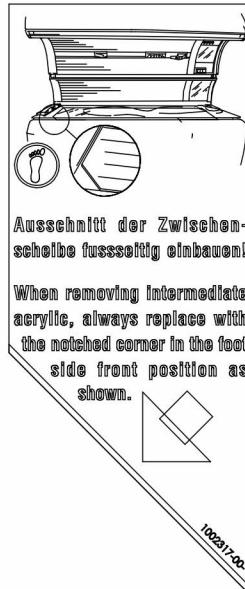

17480 / 0

11: Age warning label (1508746-..)

**ATTENTION: THIS SUNLAMP PRODUCT
SHOULD NOT BE USED ON PERSONS
UNDER THE AGE OF 18 YEARS.**

17303 / 0

The sticker is mounted on the front panel.

12: OEM parts sticker (801138-..)

ATTENTION!


Only use OEM certified replacement parts.
The use of other parts voids all
manufacturer's warranties and responsibilities.

801138-00-

13513 / 0


The sticker is mounted on the bottom of the bed base.

13: AQUA MIST sticker (801211-..)

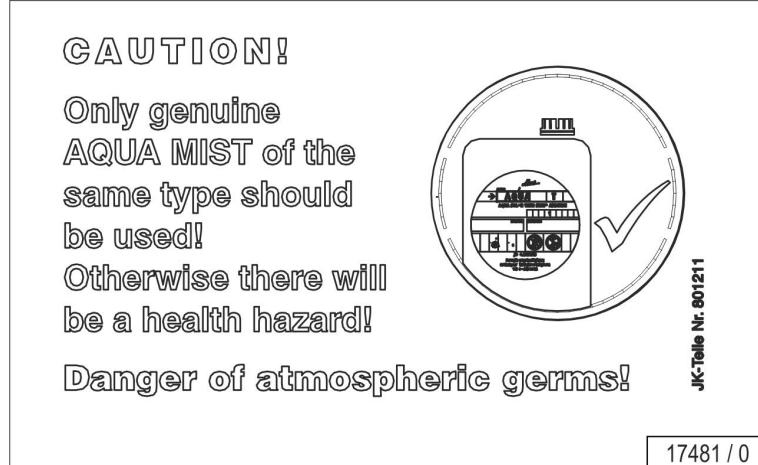

The sticker is mounted under the tray.

14: Air conditioning sticker (1002243-..)

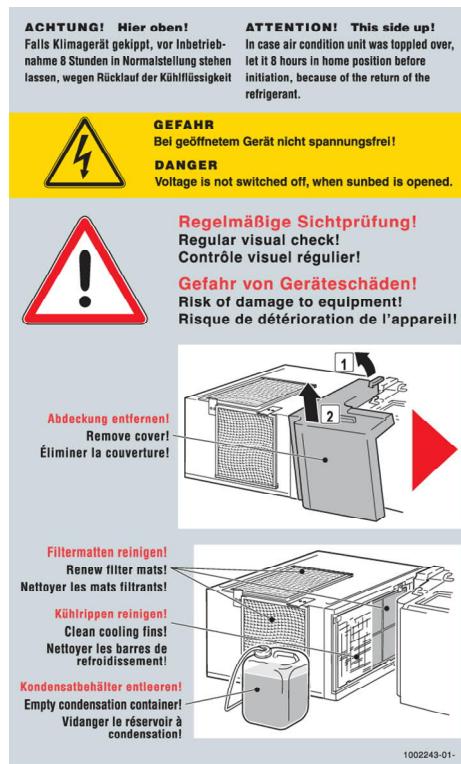

13435 / 0

The sticker is mounted on the Air conditioning.

15 : Strap sticker (1000707-..)

JK-Nr.: 1000707-..



Achtung

Erdungskabel:
Lasche nach Gerätemontage wieder montieren!

CAUTION


Ground strap:
Refit the strap after unit assembly!

17482 / 0

The sticker is attached to the ground strap.



16: Ballast support warning sticker (1003399-..)

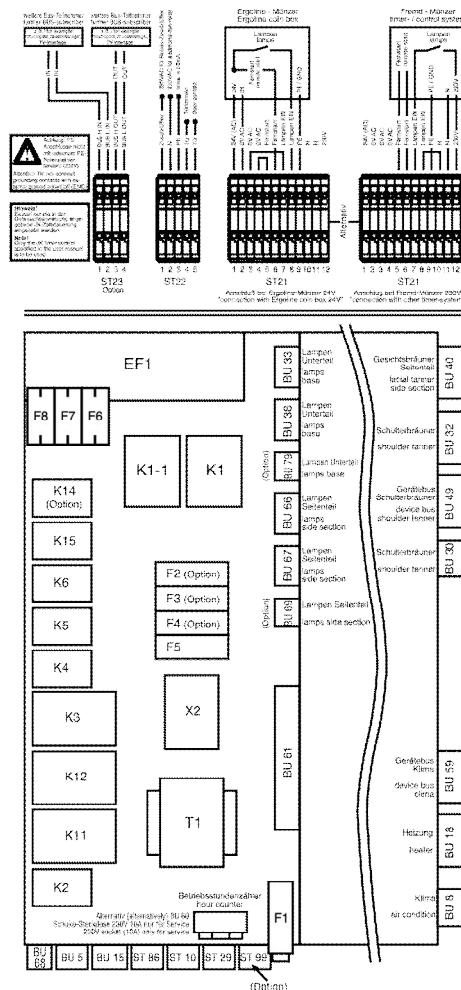

HINWEIS!
In diesem Gerät dürfen nur Originalenteile von JK-Products GmbH bzw. von JK-Products GmbH freigegebene Teile verwendet werden. Ersatzteile, die die oben genannten Voraussetzungen nicht erfüllen, jedoch den Anspruch erheben mit originalen Teilen identisch zu sein (Ident-Teile) dürfen grundsätzlich nicht eingebaut werden. Über Ausnahmen im Einzelfall entscheidet die JK-Products GmbH.

WARNHINWEIS!
Der Einbau von nicht frei gegebener Fremdelektronik oder ungenehmigte Änderungen jedweder Art an Konstruktion, Verdrahtung o. ä. führen zum Erlöschen der Gewährleistungsansprüche. Zwieiderhandlungen können **Brandgefahren** auslösen, durch die Personen infolge Frauhervergütung und/or Feuer getötet oder schwer verletzt werden können. Des Weiteren erlischt die für das Gerät ausgestellte EU-Konformitätserklärung (CE-Kennzeichen). Es besteht ferner das Gefahr, dass Sie die Ansprüche aus der gesetzlichen Produkthaftung verlieren.

NOTE!
Only genuine parts or parts approved by JK-Products GmbH or JK-Products GmbH may be used in this device. Replacement parts that do not meet the above requirements, but which, however, claim to be identical with the above-mentioned parts (pattern parts), must not be fitted under any circumstances. JK-Products GmbH will judge exceptions to this on a case by case basis.


WARNING!
The installation of unapproved electronics made by other manufacturers or unapproved changes of any kind to the design, wiring etc. shall result in any and all warranty Claims being forfeited.
Non-conformance may result in a **risk of fire**, through which, persons may be killed or badly injured as a result of smoke inhalation or fire. Moreover, the EU conformity declaration (CE certification) for the unit will be invalidated. There is also a danger of losing the right to any claims based on legal product liability.

Sicherheitshinweis:
Zum störungsfreien Gerätetrieb muß der Steuerboxdeckel immer ordnungsgemäß verschraubt sein!
Note on safety:
The control box cover must always be properly screwed on for fault-free device operation!

17483 / 0

The sticker is mounted on the power unit support.

17 : Controller connections sticker (843147-..)

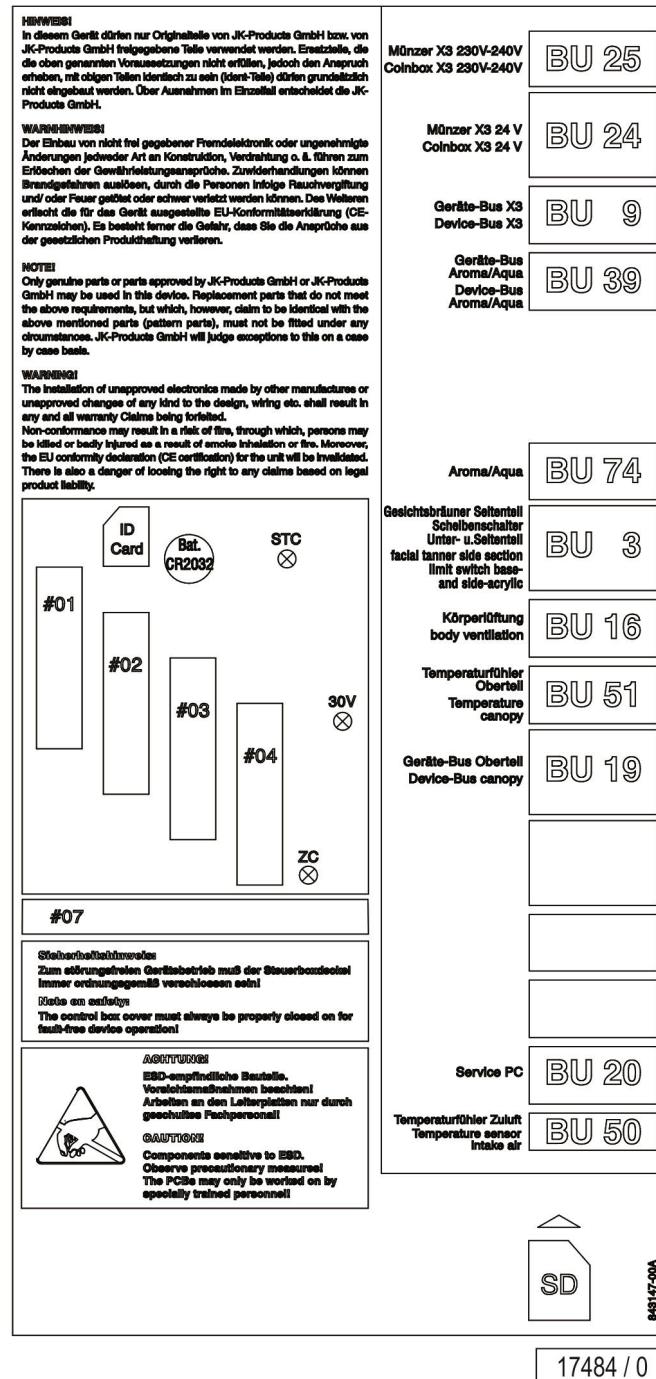

The sticker is mounted on the controller support.

18: Sticker ESD/Aquamist (1000791-..)



17485 / 0

The sticker is mounted on the AROMA canister.

19: Sticker ESD (85662-..)




17486 / 0

The sticker is mounted on several places in the device.



20: USA mains connection sticker (1511458-..)

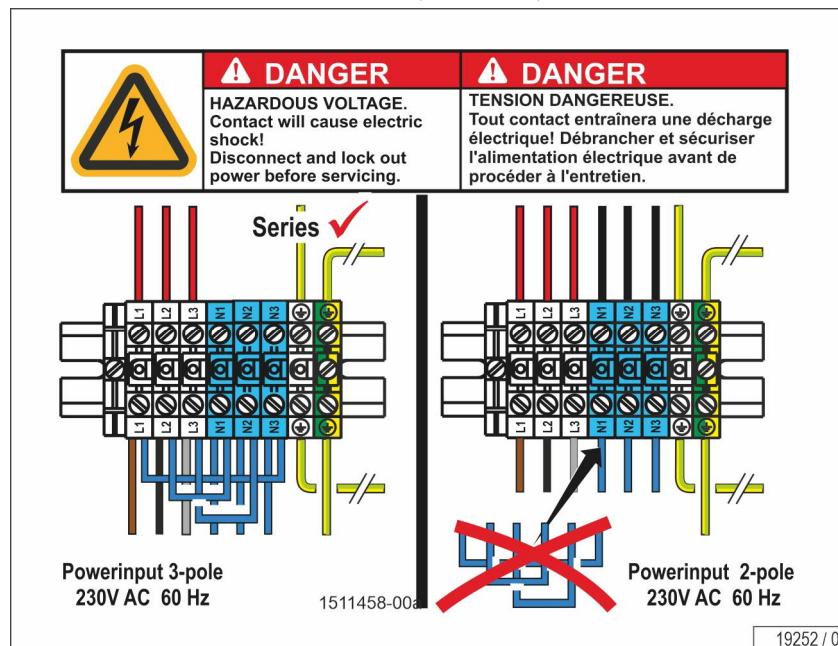

The sticker is mounted on the left side panel.

21: Canada warning sticker (801371-..)

Not applicable for this device.

22: Sticker Chip card (843148-..)

Not applicable for this device.

23: Warning sticker sensor base station (1009643-..)

Not applicable for this device.

3.10 Technical data



3.10.1 Power rating, electrical specifications and noise level

Model-No.: JK 162 / 52-4 TT HP		
Power consumption rating:	17,600 W	
Rated frequency:	60 Hz	
Rated voltage:	230 V 3 phase	230 V 2 phase ¹⁾
Rated overcurrent protection device:	70 amp 3-pole*	or 120 amp 2-pole*
Power supply:	4 x THHN AWG 6 ²⁾	or 7 x THHN AWG 6 ²⁾
Classification:	Type BF applied part lying surface (PMMA)	
Protection ingress classification:	IP20	
Protection against shock:	Class I	
Max. allowed network impedance:	Not applicable for this device.	
Rated maximum output power	normal mode: 5.26 dBm (e.i.r.p.)	
Bluetooth module:	EDR mode: 3.01 dBm (e.i.r.p.)	
Noise level 1m away from device: (without exhaust air):	62.9 dB(A)	
Noise level inside the device:	< 70 dB(A)	
Load on the acrylic panel:	max. 330 lb/150 kg	
Weight:	approx. 1,554 lb/705 kg	
Max. humidity:	< 70%	
Ambient temperature:	59-89 °F (15-32 °C)	
Max. operating altitude (above sea level):	< 6,500 ft (2,000 m)	

* Due to high inrush current caused by the electrical components used, such as ballasts, power supplies, etc., the units must be protected by slow-blow fuses.

When using circuit breakers, make sure that they have tripping characteristic D or similar.

¹⁾ For 2-phase installation, you need a special conversion box for the device to be operated. Please contact the customer service department before putting the device into operation - see page 2.

²⁾ Wire size is dependent on local electrical code. Please check with your local electrician or inspector to ensure proper size is used.



NOTE:

Adaptation of the unit is required for operating altitudes greater than 6,500 ft / 2,000 m above sea level in order to guarantee trouble-free operation.

Essential performance



NOTE:




The essential performance of the device is the UV exposure dose within specified limits, whereas the upper limit must not be exceeded. No-radiation or zero radiation does not result or contribute to an injury of the user.

Lower limit	exposure dose 0,75 MED ^{*)1}
Upper limit	exposure dose 4 MED ^{*)}

¹ *) in compliance with 21 CFR 1040.20 and FDA's timer policy



3.10.2 Dimensions

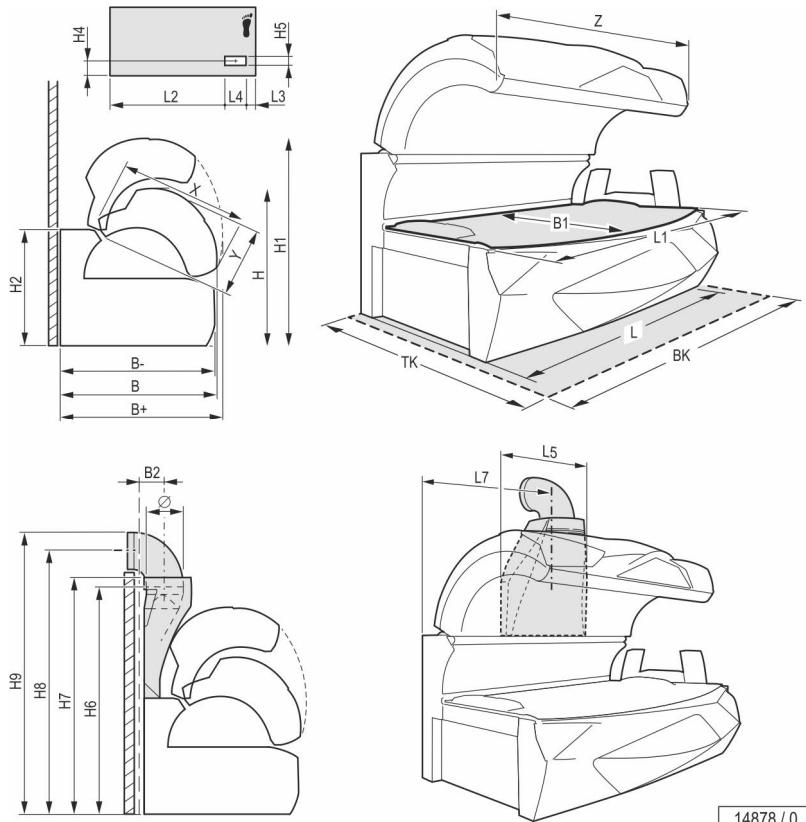

Model-No.: JK 162 / 52-4 TT HP:

B + =	60 13/16 in	1,545 mm	H =	58 11/16 in	1,491 mm
B =	59 14/16 in	1,521 mm	H1 =	76 15/16 in	1,955 mm
B - =	59 11/16 in	1,516 mm	H2 =	43 9/16 in	1,106 mm
B 1 =	37 1/2 in	953 mm	H4 =	17 11/16 in	450 mm
B 2 =	8 1/16 in	205 mm	H5 =	2 3/16 in	56 mm
L =	93 11/16 in	2,380 mm	H6 =	74 in	1,879 mm
L 1 =	86 5/8 in	2,200 mm	H7 =	78 13/16 in	2,002 mm
L 2 =	69 1/8 in	1,756 mm	H8 =	87 3/16 in	2,215 mm
L 3 =	9 1/4 in	235 mm	H9 =	92 15/16 in	2,360 mm
L 4 =	9 9/16 in	243 mm	Ø =	11 13/16 in	300 mm
L 5 =	34 1/8 in	867 mm	BK =	98 7/16 in	2,500 mm
L 7 =	44 1/8 in	1,120 mm	TK =	90 9/16 in	2,300 mm
X =	50 4/16 in	1,277 mm			
Y =	28 1/16 in	712 mm			
Z =	88 3/16 in	2,240 mm			


3.10.3 Lamp configuration

NOTICE!

The use of other lamps and filter panels voids the operating license!



Model-No.: JK 162 / 52-4 TT HP

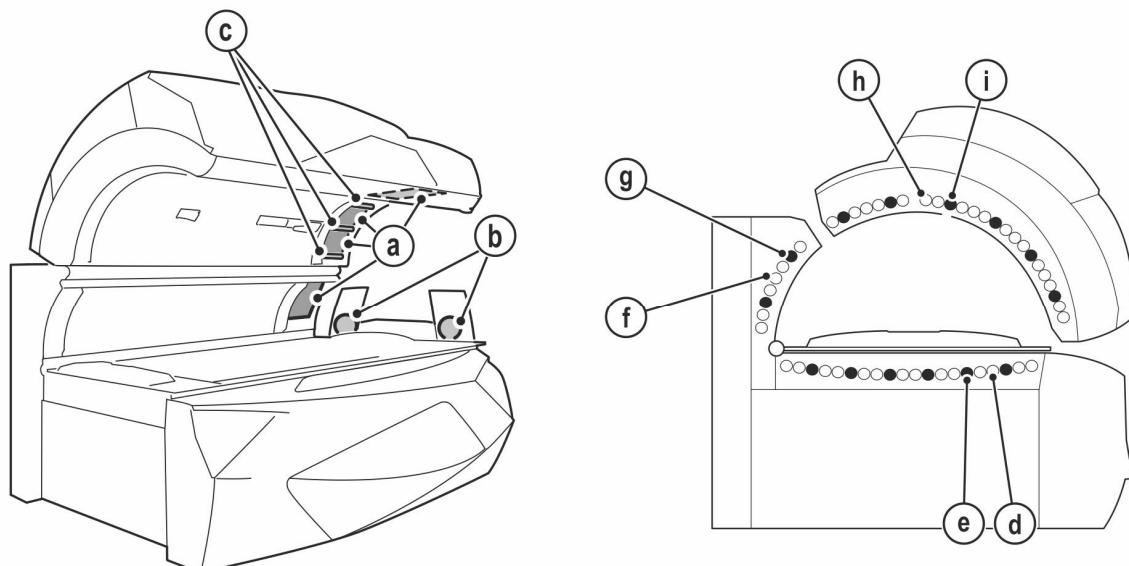

14921 / 0

Lamp set


	Number	Description	Length	Power	Control	Order-no.
a)	4	UV high pressure lamps Ergoline Ultra 520W		520 W	520 W / 420 W / 360 W / 260 W / 0 W	12645...
	4	Filter panels Ultra Performance 412H				1001526...
b)	2	UV high pressure lamps Ergoline Ultra 250W		250 W	240 W / 0 W	1500384...
	2	Filter panels Ultra Performance 912 Filter				1500635...
	2	Filter panels uncoated				1500160...
c)	3	Low pressure lamps Genesis Type R Hybrid 8W	0,31 m	8 W	8 W / 0 W	1510960...
	3	Starter S10 (5-65 W)				10047...
d)	14	UV low pressure lamps Genesis Type U HYBRID Performance 200W	1,80 m	200 W	80-200 W	1510949...
e)	6	Low pressure lamps Genesis Type R HYBRID Performance 200W	1,80 m	200 W	80-200 W	1510950...
f)	6	UV low pressure lamps Genesis Type U HYBRID Performance 200W	1,80 m	200 W	80-200 W	1510949...
g)	2	Low pressure lamps Genesis Type R HYBRID Performance 200W	1,80 m	200 W	80-200 W	1510950...
h)	18	UV low pressure lamps Genesis Type U HYBRID Performance 200W	1,80 m	200 W	80-200 W	1510949...
i)	6	Low pressure lamps Genesis Type R HYBRID Performance 200W	1,80 m	200 W	80-200 W	1510950...



3.10.4 Accent lighting

Not applicable for this device.

3.10.5 Spare parts

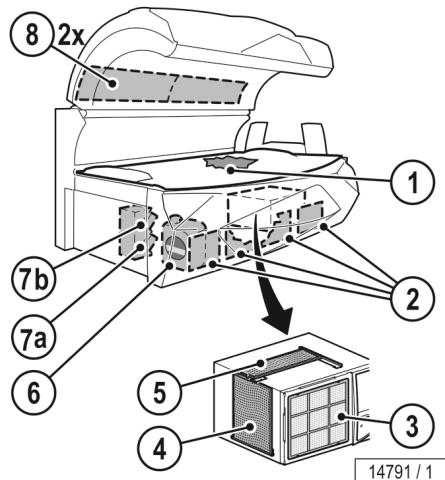

Model-No.: JK 162 / 52-4 TT HP:

Number	Description	Order-no.
1	1 Partition	1500539-..
2	4 Filter	52104-..
3	1 Filter	-
4	1 Filter	51937-..
5	1 Filter	55567-..
6	1 AQUA SYSTEM, 1.59 gal (6 l)	800617-..
7	1 AROMA, 3.38 fl oz (100 ml)	
a)	CABIN	701358-..
b)	RELAX	701357-..
8	2 Filter	843150-..
	1 Protective goggles	1027211-..
	1 AUX cable	1001855-..
	1 Allen-key	83828-..
	1 UV-lamp Type U set chip cards	1511791-..
	1 Lamp Type R set chips cards	1511792-..

For more spare parts, contact Customer Service (see page 2).

4 Why you use it



4.1 Intended use

This sunlamp product is intended exclusively for aesthetic tanning of the human skin, for one person at a time, at the age of 18 or above.

4.2 Skin condition

You may benefit from the use of this device

- if you tan in natural sunlight (skin types III, IV and V),
- if you have no acute sunburn or skin disease,
- if you take no medication or use cosmetics that may increase sensitivity to ultraviolet light.

Skin type I (Sensitive):

Invariably gets sunburn, accepts little natural sunlight and does not tan.

Skin type II (Fair):

Gets sunburn easily and severely, tans minimally or lightly and peels.

Skin type III (Average):

Often referred to as "Average" complexion, gets moderate sunburns and tans about average.

Skin type IV (Brown):

Seldom gets sunburns, tans easily and above average with each exposure.

Skin type V (Dark Brown):

This individual's skin rarely gets sunburns, tans easily and substantially.



5 Who cannot use it



5.1 Prohibited usage

5.1.1 Tanning system

WARNING!



UV radiation

Danger of skin and eye injuries or skin diseases!

The following people must not use the device:

- People under 18 years of age
- People with skin lesions or open wounds
- People with acute sunburn
- People who repeatedly suffered from serious sunburn during childhood
- People with naturally red hair and pale skin
- People with atypically discolored patches of skin
- People who have or have had skin cancer
- People who suffer from an increased risk of skin cancer (e.g. cases of skin cancer in the family)
- People among whose closest relatives a melanoma has occurred
- People who are being medically treated for photosensitivity
- People taking photosensitivity causing medicine
- People with freckles and people who tend to develop freckles
- People with more than 16 moles on their body (2 mm diameter or larger)
- People with atypical moles (atypical moles are, for example, asymmetrical moles with a diameter greater than 5 mm, with varying pigmentation and irregular edges)
- Seek medical advice if in doubt!

5.1.2 Vibra Shape (vibration plate)

Not applicable for this device.

6 What you must do to avoid serious harm



6.1 General

DANGER – HIGH VOLTAGE!



You can be killed or injured by an electric shock!

Assembly and electrical connection must comply with national and local regulations.


- The device may only be assembled, installed, upgraded or repaired by authorized, trained and qualified staff.

WARNING



Risk of electric shock!

- To avoid risk of electric shock, this equipment must only be connected to a mains supply with protective earth grounding.

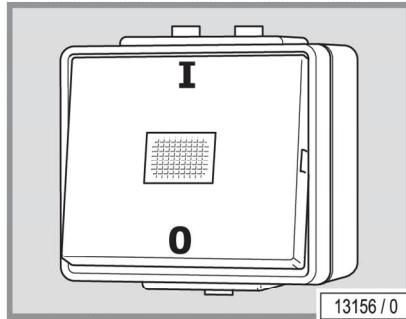

The electrical installation must have a freely accessible, all-pole disconnecting device (Master switch with clearly labelled On/Off settings, e.g. '0' and 'I' or additional indicator light) fitted externally on the site.

DANGER – HIGH VOLTAGE!



You can be killed or injured by an electric shock or injured by overexposure!

Skin damage and sunburn possible

Danger of overheating by obstruction of the air routing.

- No modification of this equipment is allowed.



6.2 Tanning system

⚠ WARNING!



Unprotected tanning may result in eye injury and eye diseases!

The unprotected eye can become inflamed on the surface, and in certain cases the retina can be damaged by excessive radiation. Cataracts can form after repeated exposure to radiation.

- Inform users to use the included protective goggles that block UV.

⚠ WARNING!



Danger of ultraviolet radiation.

Ultraviolet radiation - follow instructions.

Failure to use protective eyewear may result in severe burns or other eye injury - if discomfort develops, discontinue use and consult a physician.

- As with natural sunlight, overexposure can cause eye injury and sunburn;
- Repeated exposure may cause premature aging of the skin and skin cancer;
- Medications or cosmetics applied to the skin, may increase sensitivity to ultraviolet light;
- A person who does not tan in the sun most likely will not tan from the use of this device;
- A person having a history of skin problems should consult a physician before use;
- Overexposure should be avoided;
- The elderly, or fair skinned people who always burn easily and either never tan or tan minimally should not use this equipment.
- The appliance is not allowed to be used by more than one person at the same time.
- This sunlamp product should not be used on persons under the age of 18 years.

⚠ WARNING!



Danger of long-term consequences!

Exposure to ultraviolet radiation may increase the likelihood of developing skin cancer and cause serious eye injury.

- Overexposure should be avoided.
- This sunlamp product should not be used on persons under the age of 18 years.

WARNING!



WARNING!

This device emits electromagnetic fields!

EM-Emission inside device:

- 5 Hz – 2 kHz max. 2.1 μ T (1.5 kV/m)
- 2 kHz – 400 kHz max. 3.3 μ T (1.3 kV/m)
- 9 kHz – 30 MHz max. 35 db μ V/m (56 μ V/m)
- 30 MHz – 1 GHz max. 95 db μ V/m (56 μ V/m)
- 1 GHz – 6 GHz max. 90 db μ V/m (32 μ V/m)



If you have an active implant or body worn medical device (e.g. insulin body pump, blood glucose sensor) consult your health care provider before use of this device.

- Seek medical advice if in doubt!

6.3 Sensor

Not applicable for this device.

7 What you must do to avoid other harm



7.1 General

WARNING!



WARNING!

Infections can be transferred through skin contact.

All items/equipment components, with which the user can come into contact during tanning, must be cleaned and disinfected after every tanning session:

- Protective goggles
- Acrylic panels
- Handles and control panel
- AUX cable
- Adjustable air nozzles
- Use only EPA registered quaternary ammonium germicidal based detergent solution (follow product label for use-dilution), such as Lucasol One Step, to clean and disinfect the tanning system and its components.

7.2 Tanning system

WARNING!



WARNING!

Danger of long-term consequences!

- Overexposure should be avoided.
- The elderly, or fair skinned people who always burn easily and either never tan or tan minimally should not use this equipment.
- This sunlamp product should not be used on persons under the age of 18 years.
- While there is no immediate clinical evidence of UV exposure and its effects upon expectant mothers, it is strongly advised that expectant mothers be discouraged from using the tanning unit.

7.3 Sensor

Not applicable for this device.

7.4 Vibra Shape (vibration plate)

Not applicable for this device.



7.5 Foreseeable misuse

The following uses of the device are explicitly prohibited:

- Only one person at a time may use the device.
- People with body weight of more than 330 lb (150 kg) may not use the device.
- Exceeding the specified tanning times may result in burns and long-term skin damage.
- When the device is in use no other people may remain in the room. This includes both children and adults.

7.6 Obligations of the studio operator

- As the studio operator you are responsible for informing users on possible risks of tanning.
- All danger and safety labels attached to the device must be observed! Do not remove them!
- No safety equipment (e.g. panel switch) may be removed, deactivated, or otherwise bypassed as it can impair the safe operation of the device!
- The device may only be used when in completely operational condition!
- It is the studio operator's responsibility to ensure the tanning system is handled properly.

7.7 Staff training

As the operator you are obligated to train your personnel and instruct them in the established legal and accident prevention regulations. Ensure that your staff have understood and observe these instructions. This includes:

- information regarding intended use and foreseeable misuse
- safety instructions
- operating instructions
- cleaning/disinfection and maintenance instructions
- Please make the instructions and additional information from the manufacturer available to the staff at their workplace.
- Regularly check your staff's safety consciousness and awareness of the risks, taking into account your operating instructions.
- Please ensure that the cleaning/disinfection intervals are complied with, in particular cleaning and disinfection after every use.
- Please ensure that your staff provide your customers with appropriate advice, in particular regarding tanning times, possible risks and the operation of the device.

8 Risks of using this device



⚠️ WARNING!



Danger through ultraviolet radiation!
Danger of skin and eye injuries or skin diseases!

- This sunlamp product must not be used on individuals who have had skin cancer or have a family history of skin cancer.
- Persons repeatedly exposed to ultraviolet sunlamp products must be evaluated for skin cancer regularly.
- This sunlamp product must not be used by people among whose closest relatives a melanoma has occurred.

Contraindications! The following people must not use this system:

- People under 18 years of age
- People with skin lesions or open wounds
- People with acute sunburn
- People who repeatedly suffered from serious sunburn during childhood
- People with naturally red hair and pale skin
- People with atypically discolored patches of skin
- People who have or have had skin cancer
- People who are being medically treated for photosensitivity
- People taking photosensitivity causing medicine
- People with freckles and people who tend to develop freckles
- People with more than 16 moles on their body (2 mm diameter or larger)
- People with atypical moles (atypical moles are, for example, asymmetrical moles with a diameter greater than 5 mm, with varying pigmentation and irregular edges)

9 Benefits of using this device



Traditional UV light with the addition of (beauty) red light results in a natural looking tan. Optimal tanning results are achieved due to the synergistic combination of light spectrums.

10 How to decide to use this device



This device is intended only for commercial use and not for private use.

Some people are not permitted to use the device, see chapter 5.

The device may not be used by people with reduced physical, sensory or mental capabilities or who have insufficient experience and/or knowledge without supervision. A person responsible for their safety must ensure, either by supervision or instruction, that the device is used correctly and safely. If there is any doubt, use of the device by such categories of persons is to be forbidden!

11 What to do before the treatment



Please check for any conditions that may prohibit the use of this device – see chapters 5, 6, 8, 10.

Make-up in the tanning device

Please do not wear make-up. Clean skin absorbs UV light better.

A wide variety of materials is found in cosmetics. Whether these are emulsifiers, greases, or fragrances - in conjunction with UV radiation, these can cause skin allergy reactions.

Remove all make-up. This is essential. Make-up clogs the pores in the face. The UV light from a tanning device re-opens these clogged pores so that not only does the light penetrate the skin, but also the allergenic components of the make-up.

A further negative consequence of make-up is that the appearance of the skin deteriorates irreversibly over a period of time. However good make-up looks - in conjunction with UV light, it causes more harm than good. Therefore remove all make-up before a session in your tanning device - it will be much better applied afterwards.

- Remove all jewelry (including piercings) prior to the application.
- Remove cosmetics before using the device and do not use sun screen products.

Are medications and the tanning device compatible?

It is known that some medications increase the skin's sensitivity to UV light. This probability is particularly high with antibiotics, sulphonamides, psychiatric drugs, tranquillizers, anti-diabetes treatments and diuretics. Sun creams which contain psorale or kumarin also increase skin sensitivity. If you are in any doubt, first ask your doctor.

Contact lenses in the tanning device?

Like any other tanning device user, contact lens wearers should also wear special protective goggles which protect the eyes against UV radiation.

Even when using contact lenses with UV protection, the protective goggles are mandatory.

Showering after leaving the tanning device?

The browning is in the skin, not on it. Therefore you cannot wash it off. Treat your skin to a moisturizing cream after showering.

12 How to use this device



12.1 How to set this device up

After commissioning and cleaning/disinfection, the device is ready for operation.

Via the Control Center the most important settings can be made before starting the device. During the tanning session, functions may be activated or deactivated via the control panel inside the device.

See chapter 12.5.

12.2 How to make sure the device works

Before using the device, check for damage to filter panels and acrylic panels.

12.3 How to use it to treat your condition

- No other persons are permitted inside the cabin while the device is in use.
- UV light can cause irreversible skin or eye damage. The unprotected eye can become inflamed on the surface, and in certain cases the retina can be damaged by excessive radiation. Cataracts can form after numerous repeated tanning sessions. Use the UV-impermeable protective goggles provided (order no. 1027211-..). Contact lenses and sun glasses are no substitutes for protective glasses.
- Always comply with the recommended application times.

12.4 How to choose correct treatment regimen

12.4.1 Skin types

Skin type I (Sensitive):

Invariably gets sunburn, accepts little natural sunlight and does not tan.

Skin type II (Fair):

Gets sunburn easily and severely, tans minimally or lightly and peels.

Skin type III (Average):

Often referred to as "Average" complexion, gets moderate sunburns and tans about average.

Skin type IV (Brown):

Seldom gets sunburns, tans easily and above average with each exposure.

Skin type V (Dark Brown):

This individual's skin rarely gets sunburns, tans easily and substantially.



12.4.2

Tanning times



NOTE:

The tanning times given in the table must not be exceeded!



DANGER ULTRAVIOLET RADIATION						
Follow all instructions. Avoid over-exposure as with natural sunlight. Over-exposure can cause eye and skin injury and allergic reactions. Repeated exposure may cause premature aging of the skin and skin cancer.						
MAXIMUM TANNING TIME: 10 MINUTES						

Recommended Exposure Schedule

Skin Type	Week 1 1st - 3rd Treatments	Week 2 4th - 6th Treatments	Week 3 7th - 10th Treatments	Week 4 11th - 15th Treatments	Weekly Subsequent Treatments
I - Sensitive	Not Advised	Not Advised	Not Advised	Not Advised	Not Advised
II - Fair	2 Min.	4 Min.	6 Min.	8 Min.	10 Min.
III - Average	2 Min.	4 Min.	6 Min.	8 Min.	10 Min.
IV - Brown	2 Min.	6 Min.	8 Min.	10 Min.	10 Min.
V - Dark Brown	2 Min.	6 Min.	8 Min.	10 Min.	10 Min.

WEAR PROTECTIVE EYEWEAR! Failure to do so may result in severe burns or long term injury to the eyes.

Certain medications or cosmetics may increase your sensitivity to the ultraviolet radiation. Consult a physician before using sunlamps if you are using medication, have a history of skin problems or believe yourself to be especially sensitive to sunlight. If you do not tan in the sun, you are unlikely to tan from the use of this product.

This unit utilizes UV lamps. Replace only with Ergoline ULTRA 520W (high-pressure / facial tanner) and Ergoline ULTRA 250 W (high-pressure / shoulder tanner) and Genesis Type U HYBRID Performance 200W or other recognized equivalents.

To obtain the recommended exposure distance, lie in the middle of the acrylic glass panel and lower the canopy to the closed position. Canopy will stop at the correct distance from lounge acrylic surface.

Warning: The use of any other position may result in overexposure!

Tanning can begin on a regular basis. A tan normally appears after several exposures and is maximized after four (4) weeks of exposure in accordance with the recommended exposure schedule for your skin type. In order to maintain your tan 10 minute sessions are recommended from the fifth week onward.

Do not tan more than once in a 24-hour period.

Use protective eyewear whenever unit is energized. Read and follow instruction booklet before using this product. To avoid electrical shock during maintenance disconnect power from unit. This product is in conformity with performance standards for sunlamp products under

21CFR1040.20, FDA.


Part No: 1019217-XX

18534 / 2

The sticker is mounted on the lamp cover foot end.

12.5

How to operate the device

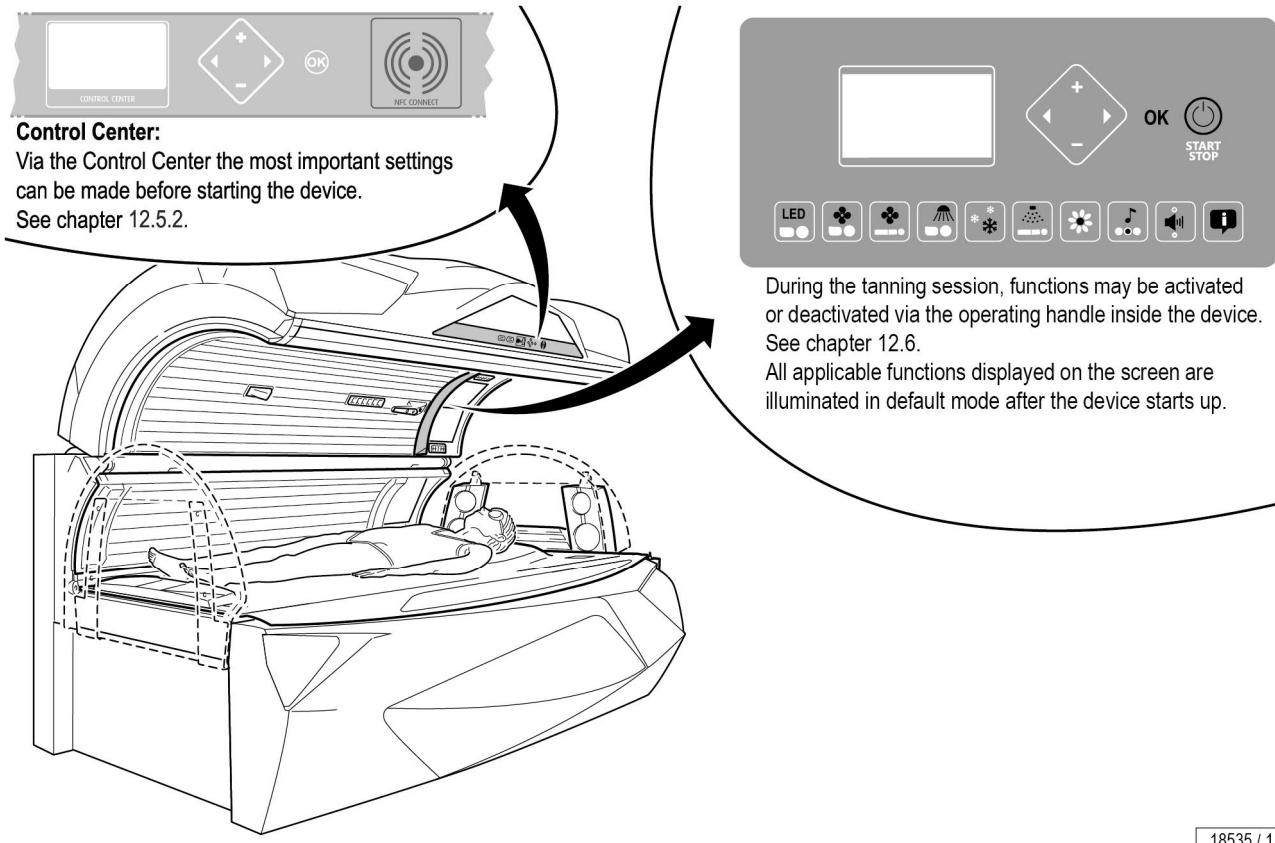


During the tanning session, functions may be activated or deactivated via the operating handle inside the device. See chapter 12.6.

All applicable functions displayed on the screen are illuminated in default mode after the device starts up.

18535 / 1

12.5.1

Control Center

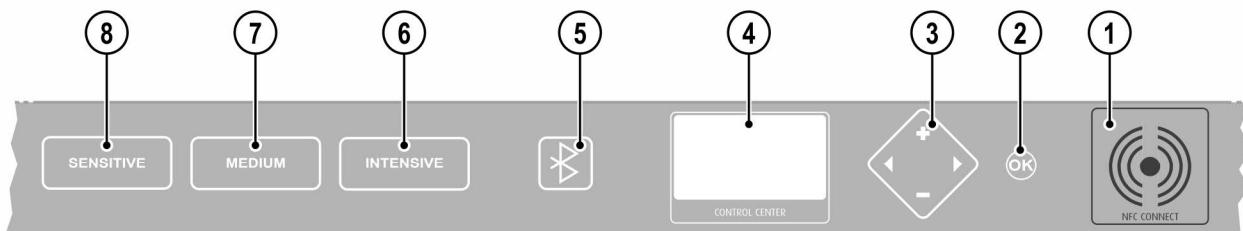

17977 / 0

1. NFC Connect
2. OK button (not always visible)
3. Feature selection
4. Display – Functions ...
5. Bluetooth® functions (only before tanning)
6. Application type: INTENSIVE
7. Application type: MEDIUM
8. Application type: SENSITIVE



12.5.2 User settings overview

User settings overview

All of a device's available functions can be pre-set on the control panel and Control Center. Fine adjustment is carried out during use.



NOTE:

Pressing and holding for 3 seconds opens the 'User settings' overview.

If the START/STOP button is pressed and held for 3 seconds while in standby mode, accent lighting and music switch off or back on.

Function	Pictogram selection / key sequence	Description / display
Body fan	HIGH AUTO LOW	Pre-select low or high cooling or set to automatic.
Face fan	HIGH AUTO LOW	Pre-select low or high cooling or set to automatic.
Facial tanner	MAX MED OFF	Select maximum or medium power or switch off.
Hybrid Performance LED	LED 3	Setting the intensity or switching off the Hybrid Performance LEDs.
Tempronic/Air conditioning	WARM AUTO COOL	Pre-select low or high temperature or set to automatic.
AROMA	 	Switch AROMA on or off
AQUA SYSTEM	 	Switch the AQUA SYSTEM body and head nozzles on or off together or individually.
Shoulder tanner	 	Switch the shoulder tanner on or off
Voice Guide	 	Switch the Voice Guide on or off
Volume	HIGH AUTO LOW	Basic volume setting: Loud, quiet or automatic.
Audio system	AUX SALON SOLARIUM OFF AUTO	Select the audio system or switch off
Welcome	WELCOME	
	Welcome text in standby mode	




NOTE:

If a button is not pressed within 30 seconds, the settings return to default mode.

The following pre-settings can be selected as standard ex works in the Control Center and control element. Functions can also be activated in the configuration menu (see Professional Setup Manager 1014640-..).


	Standard	Additional	Note
	HIGH ► AUTO LOW	yes	-
	HIGH ► AUTO LOW	yes	-
	✓ ► ✗	yes	-
			Settings by the user
	►	yes	-
	► ✓ ✗	yes	-
	WARM ► AUTO COOL	no	ja
	► ON OFF	no	yes
	HIGH ► AUTO LOW	no	yes
	► PERSONAL SALON SOLARIUM OFF AUTO	no	yes
	3	no	yes
WELCOME XXXXXX	no	yes	Input 16 characters (e.g. name)

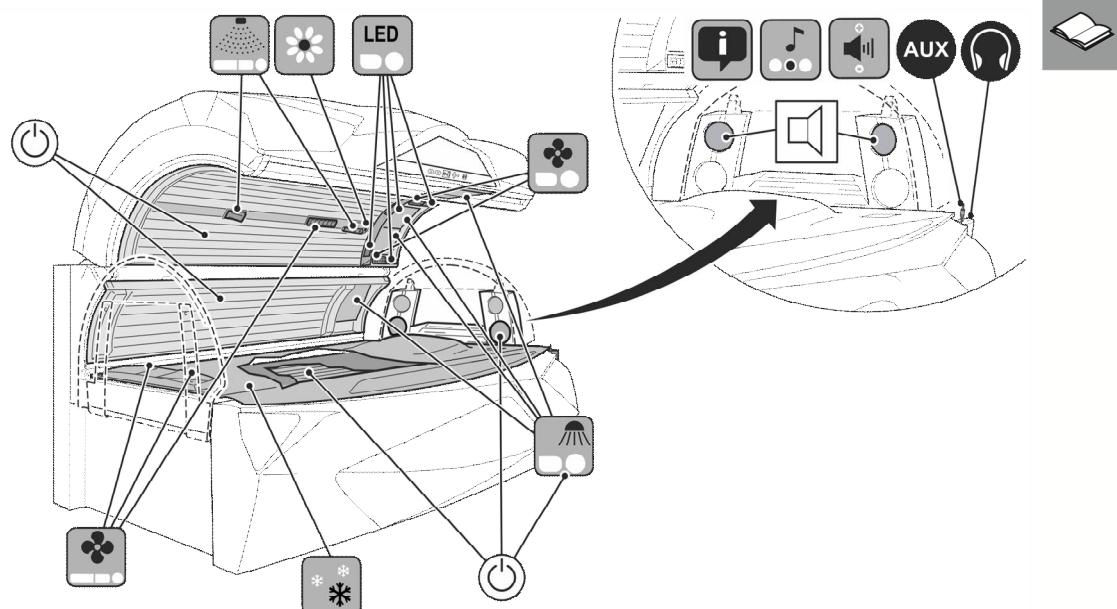

17975 / 1



Control panel

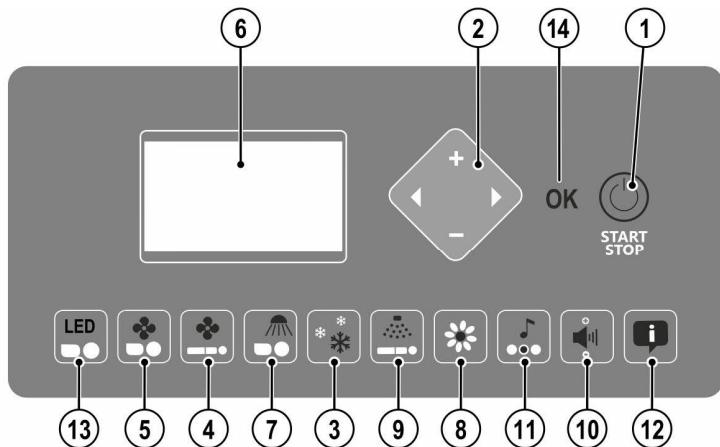

17976 / 0

1. START/STOP button
2. Feature selection
3. Temptronic pictogram
4. Body fan pictogram
5. Facial fan pictogram
6. Display – Functions ...
7. Facial tanner pictogram
8. AROMA SYSTEM pictogram
9. AQUA SYSTEM pictogram
10. Volume pictogram
11. Audio system/channel pictogram
12. Voice Guide pictogram
13. Hybrid Performance LED pictogram
14. OK button (not always visible)



NOTE:

All applicable functions displayed on the screen are illuminated in default mode after the device starts up.




12.6.1 Navigation

Via the Control Center and control panel the most important settings can be made before starting the device. During the tanning session, functions may be activated or deactivated via the control panel inside the device.

Navigation unit - Control Center

The Control Center enables comfortable access to special functions.

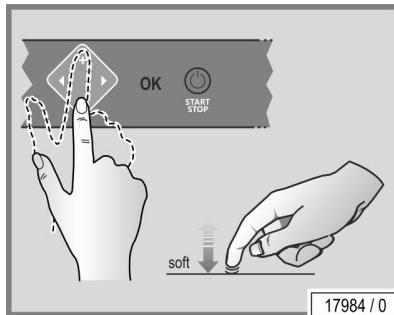

Acknowledge selection.

This button's availability is displayed on the Control Center, depending on the function that has been selected.



Select device function.



Functions or settings can be activated, deactivated or changed.



NFC¹-Point

Data exchange via contactless near field communication – individual settings can be stored on a NFC medium (e.g. NFC tag) and called up again during the next session.

Information about the setting options can be found in the separate user manual 'Control Center' (order-no.: 1014640-..) in the folder Technical Documentation.



NOTE:

Only for the studio operator:

To open the Professional Setup Manager, press button  for approximately 3 seconds in standby mode.

All of the pictograms displayed on the screen are illuminated in default mode after the device starts up.

The relevant pictogram will start to flash after a device function has been selected. The settings or selection options available for the selected function will be displayed on the screen.



NOTE:

Pressing and holding  button for 3 seconds activates the demo mode. This allows you to demonstrate the following functions without tanning.






NOTE:

If no button is pressed within 8 seconds, the settings return to default mode.


¹ Communication system: NFC (Near Field Communication); ISO/IEC 14443A/MIFARE max. operating range: approx. 50 mm, frequency range: 13.56 MHz


12.6.2 Description

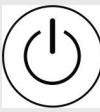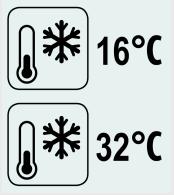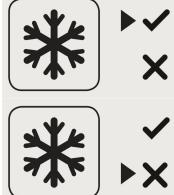Pictogram selection / Navigation key sequence	Description / Display	Information
 	START/STOP Switches the device on / off during the tanning session The process is interrupted.	The tanning time continues to run if the device is switched off during tanning session.
 	Function selection face fan, body fan, facial tanner	
 	Regulate the temperature, -16°C (~61°F) - ~32°C (~90°F) (Climatronic Plus) -18°C (~64°F) - 30°C (86,0°F) (Climatronic)	Optional
 	Air conditioning failure (Climatronic) To switch the air conditioning unit (Climatronic) off: Press the minus button for 3 seconds To switch the air conditioning unit on: Press the plus button for 1 seconds	Optional E005: If the automatic air conditioning (Temptronic) fails, the acrylic panel is cooled without adjustment control. The user can only switch the cooling system on or off in this case. During tanning the air conditioning can only be switched back on after 3 minutes.
 9  1		Regulate the body fan, levels 1-9
 8  0		Regulate the face fan, levels 0-8
 3  2  1  0		Regulate the facial tanner: levels 0-3 Switch the facial tanner off: Press the minus button for 3 seconds Switch the facial tanner on: Press the plus button for 1 seconds




The face cooling air nozzles can be adjusted manually.

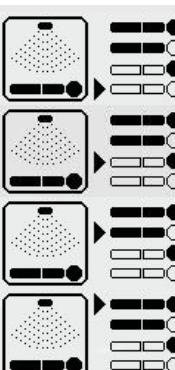Pictogram selection / Navigation key sequence	Description / Display	Information
   	<p>Switch on AROMA: Press the plus button for 1 second.</p> <p>Switch off AROMA: Press the minus button for 1 second.</p>	
 	<p>AQUA SYSTEM</p> <p>Switch the AQUA SYSTEM head nozzle on: Press the plus button.</p> <p>Switch the AQUA SYSTEM body nozzle on: Press the plus button.</p> <p>Switch the AQUA SYSTEM off: Press the minus button.</p>	
 14 	<p>Adjust volume: levels 0 - 14</p>	
 10		
 1		
 0		
  	<p>Audio system</p> <p>SD card selection: Press the plus button</p>	<p>The +/- buttons can be used to select titles for internal music in demo mode and during tanning.</p>
 	<p>Studio channel selection: Press the plus button</p>	
 	<p>Select own music (MP3 player)</p> <p>AUX or Bluetooth®: Press the plus button</p>	
 	<p>Switch off the audio system: Press the minus button</p>	<p>The +/- buttons can be used to select titles for internal music in Bluetooth® mode and during tanning.</p>

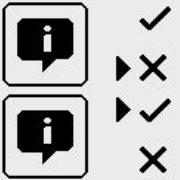Pictogram selection / Navigation key sequence	Description / Display	Information
 	 Switch off Voice Guide: Press the minus button	The voice guide provides acoustic device operation assistance and is only available if an audio system is fitted. Voice Guide is switched on by default during stand-by and tanning.
 	 Regulate the Hybrid Performance LED, levels 0-3	

12.7 Selection of the application

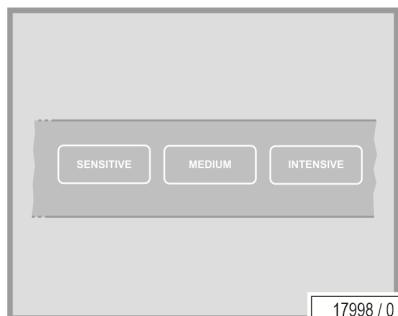

Before the session begins, the desired application type must be selected. There are 3 possible selections:

- **INTENSIVE:**
100% power of Genesis Type U HYBRID Performance¹ and Genesis Type R HYBRID Performance² and facial tanner
- **MEDIUM:**
approx. 85% power of Genesis Type U HYBRID Performance and Genesis Type R HYBRID Performance, facial tanner level 2
- **SENSITIVE:**
approx. 65% power Genesis Type U HYBRID Performance and Genesis Type R HYBRID Performance, facial tanner level 1
- Press the button for the desired type of application.



NOTE:

In each program, the Hybrid Performance LEDs start with 100% performance.

During the tanning session, the intensity of the facial tanner and the Hybrid Performance LEDs can be changed at any time, or switched off completely. If no selection of the application type is carried out, the device starts up with the application type „INTENSIVE“

¹ UV low pressure lamps

² Low pressure lamps



12.8 Operating display of the sensor

Not applicable for this device.

12.9 Sun Select setting with the outer interface

Not applicable for this device.

12.10 Measurement (automatic mode)

Not applicable for this device.



12.11 Connections

⚠ WARNING!



WARNING!

Danger of hearing damage!

Using earphones or headphones at high volume can result in permanent hearing damage.

- Do not set the volume too high.

12.11.1 Couple with Bluetooth® devices

The user can connect their own private MP3 player or smartphone with Bluetooth® to the device.

Coupling between the solarium and the private smartphone or MP3 player of the user is possible via Bluetooth®¹ Connect or with the stereo jack plug. When using Bluetooth Connect, the user can transfer their favourite music into the solarium wireless (e.g.) from their smartphone. The operation is carried out during the tanning session conveniently via the Control Center.



NOTE:

- Ergoline cannot guarantee full compatibility with all Bluetooth® devices.
- Before coupling a device with this system, read the user manual to find out more about the Bluetooth® compatibility.
- Make sure that the Bluetooth® function is activated on your device and that the device is visible for all other Bluetooth® devices.

Connect to a Bluetooth® MP3 player

Activate the Bluetooth® function and search for Bluetooth® devices.

Activate the Bluetooth® function on the iPhone®:

iPhone®:

Settings → General → Bluetooth® → ON (search for units is started.)

Activate the Bluetooth® function on the Android® smartphone:

Android® smartphone:

Press the Bluetooth® symbol → SCAN (search for units is started.)

¹ Communication system: IEEE 802.15 WPAN Bluetooth® class II version 3.0 compatible with Bluetooth version 2.1+EDR, 1.2, and 1.1; max. operating range: approx. 10 m; frequency range: 2,4 GHz range (2,4000 GHz – 2,48 GHz); compatible Bluetooth® profiles: A2DP, AVRCP; equipment supported at the same time (Pairing): one device




NOTE:

The tanning device can store a maximum of ten coupled Bluetooth® devices. In the event of recoupling the tanning device must first be deleted from the Bluetooth® device. The activation process must then be repeated.

Activate Bluetooth® Connect on the Control Center



- Press the Bluetooth® symbol on the operating panel.

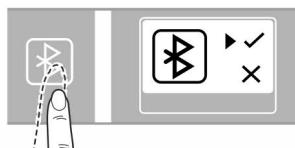

- Activate Bluetooth® Connect on your audio device.
- Select e.g. ERGOLINE 01.




Successful coupling on the operating panel.



NOTE:

Start with the audio player music on the Bluetooth® MP3 player.

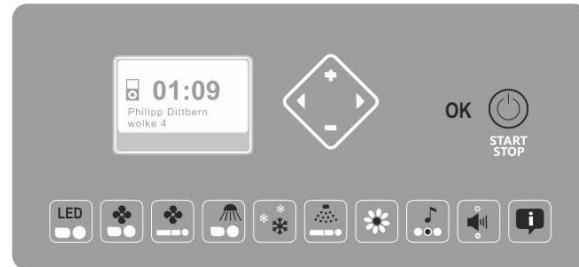

17982 / 0

Display in the control panel in a successful coupling.



NOTE:

When successfully coupled, the Bluetooth® symbol lights up continuously blue.

12.11.2 MP3 player and smartphone

⚠ CAUTION!



The heat might damage your MP3 player or smartphone!




- Do not place the unit directly in the radiation area.

The cable might tear, the connectors might be damaged.

- Do not pull on the cable when unplugging the MP3 player or smartphone.

The user can connect their own private MP3 player or smartphone with Bluetooth® to the device.

1. Jack plug (JK Part # 4630911001-.) provided for use by the studio operator (2 x 3.5 mm stereo jack plug).
2. MP3 player or smartphone connection: see the documentation from the manufacturer of the MP3 player or smartphone.

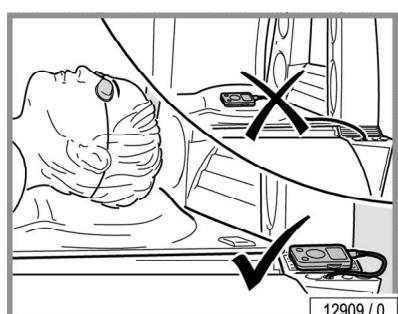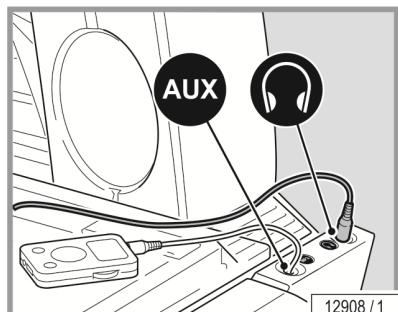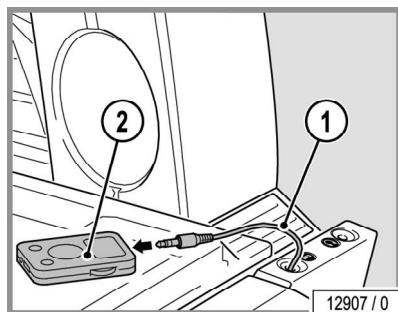

Connecting:

- Insert the AUX cable (stereo jack plug) in the audio system.
- Plug the AUX cable into the MP3 player or smartphone.
- Select the music title and then start.
- Where necessary, set the volume level on the MP3 player or smartphone so that the sound signal from the MP3 player or smartphone is recognized by the audio system.
- Pre-setting: approx. 75% of the volume on the MP3 player
- Coupling a smartphone with the device via Bluetooth®: see page 58

NOTE:

The personal use of an MP3 player during tanning is not a public broadcast in copyright terms, so that there is no obligation for registration or payment by the studio owner to the corresponding copyright owning company. Copyright laws apply to the owner of the MP3 player.

The same regulations apply for the official use of MP3 music files as for other music sources: As studio owner you are only allowed to use original audio obtained from legal sources. Proof of the legitimacy of any MP3 audio used should be available for viewing upon demand of any regulatory body. Please refer to <http://www.riaa.com> for more information.



12.12 How to use



NOTE:

Before use the device must always be cleaned and disinfected! See chapter 12.15 How to clean/disinfect.

12.12.1 How to start

The device is started automatically with an internal or external timer.




NOTE:

See separate document 'Timer Control Module' #1008370-.. in the Technical Documentation folder.

Preparation time

Before the lamps turn on, you have a few minutes to prepare yourself. During the preparation time you can switch on the lamps at any time by pressing the START/STOP button.

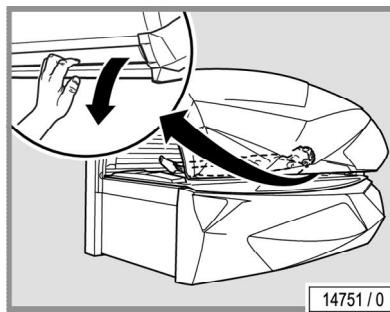

14751 / 0

- Remove your clothes and all jewelry.
- Remove make-up and other cosmetics.
- Remove glasses. Contact lenses can be worn under the protective goggles.
- Put on the protective goggles.
- Lie down on the bed surface.
- Close the canopy.
- Press the START/STOP button.

12.12.2 Application tips

Tanning device

When describing the operation, the tanning device is often shown with the canopy open. During tanning the canopy must always be closed.

Proper operation is not possible with the canopy open.

- Lie down and place your head on the headrest or as indicated by the ergonomic base panel to maintain the proper distance between skin and lamps.

12.12.3 How to stop

The tanning session ends automatically after the preset time has run down. The internal body cooling fan continues to run for approx. three (3) minutes to sufficiently cool the device.

- To interrupt a tanning session, press and hold button START/STOP any time during a session.
- Press START/STOP again to restart the lamps.



12.13 What to do after the treatment

12.13.1 Cleaning and disinfection

WARNING!



Infections can be transferred through skin contact. All items/equipment components, with which the user can come into contact during tanning, must be cleaned and disinfected after every tanning session:

- Protective goggles
- Acrylic panels
- Handles and control panel
- AUX cable
- Adjustable air nozzles
- Use only EPA registered quaternary ammonium germicidal based detergent solution (follow product label for use-dilution), such as Lucasol One Step, to clean and disinfect the tanning system and its components.



NOTE:

For thorough cleaning and disinfection, the product must be allowed to remain in contact with the surface being cleaned and disinfected for a specified amount of time. This specified amount of time depends on the cleaner/disinfectant which is used. Follow the directions for use of the cleaner/disinfectant.

After this specified application time remove remaining cleaner/disinfectant from the surfaces of the device. Use a single-pass paper towel.

Once the remaining cleaner/disinfectant has been removed, perform a visual inspection for signs of any remaining cleaning solution, or areas still soiled. If either are remaining, repeat the cleaning process until a visual inspection can be successfully completed without the need for additional cleaning.

If any of the components listed above cannot be successfully cleaned and disinfected, they must be replaced.





12.14 Safety instructions for cleaning/disinfection and maintenance

DANGER – HIGH VOLTAGE!



You can be killed or injured by an electric shock!
During maintenance and service work there is a risk of touching live parts and receiving a fatal electric shock.

- The device must not be in use.
- Disconnect device from power supply before opening.
- Do not reach behind covers and do not open the device with the device switched on.
- Refit any safety devices (e.g. switches) after the work has been completed.

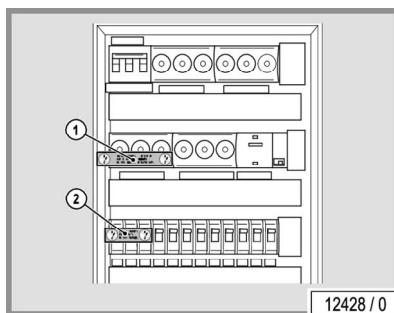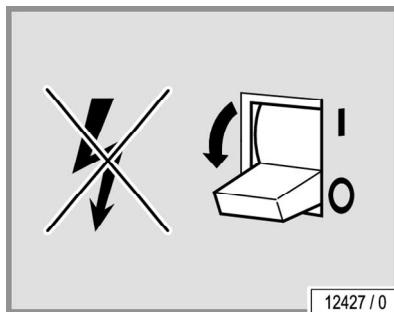

If work is to be performed on the devices, they must be deenergized. That means that all live cables must be switched off.

Switching off the device alone is insufficient, as certain points may remain energized. Therefore, switch off fuses before performing any work and – if possible – remove them.

Accidental reactivation can result in serious accidents. Immediately after deenergizing, secure all switches or fuses used for deenergizing against being switched on again.

- Lock the fuse box with a padlock.

In case of circuit breakers which cannot be unscrewed, a strip of adhesive tape may be applied over the actuation lever (1 + 2) with the wording 'Do not switch, danger!'.
A prohibition sign with the text:

'Work in progress!'

'Place:'

'Sign may only be removed by:'

must always be securely mounted immediately.

It is imperative that you, the owner and/or operator, check all components for any signs of wear which may occur through normal use.

- Observe the cleaning/disinfection and maintenance schedule specified in the tables in chapter 12.15.2 and 12.16.2.
- Only replace low pressure lamps with the original type or compatible lamps.

Should you require any assistance or need additional information, please call Global Service USA at 800 445 0624, Monday through Friday 8:00 a.m. - 5:00 p.m. CT.



NOTE:

The service of the device may only be carried out by a (Level 2) JK Factory Certified Technician.



⚠ WARNING!



Injuries and damage due to improper maintenance!

If incompatible UV lamps are installed, the share of harmful UV radiation may increase. Sunburn and other skin damage may result.

Improper maintenance can cause the device to be damaged or to become unusable.

- Always comply with the specified maintenance intervals, see chapter 12.16.2 Maintenance plan.
- Follow the cleaning/disinfection instructions.

Danger of burns from hot lamps and device parts!

If the tanning device was operated before being switched off and opened, lamps and housing parts may be hot.

- Before replacing lamps or other spare parts, wait until lamps and metal parts of the device have cooled down sufficiently.

Danger of squeezing!

When opening and closing the panels and acrylics, you can squeeze your fingers.

Maintenance and care play a major role in deciding the devices performance and appearance over its useful life. Diligent compliance at the specified intervals is absolutely necessary to ensure the proper operation of your device.

The current operating hours of the parts requiring maintenance can be displayed in the Presetting mode – see separate document 'Professional Setup Manager' #1014640-.. in the Technical Documentation folder.

Please note that the device must be inspected regularly by a JK authorized representative every 12 months (from date of installation) to maintain its proper condition!

Acrylic panel

NOTICE



Cosmetics may cause cracks to form on the surface!

Over the long-term cosmetics or sunscreens cause damage (e.g. crack formation on the surface).

- Ensure that cosmetics and sunscreens are removed before tanning.

The acrylic panels for the tanning devices are produced of acrylic glass developed especially for this application. The acrylics used feature a particularly high UV permeability and resistance, as well as an easy-care, hygienic surface that is gentle on the skin.

The acrylic panels are formed to device-specific shapes in a technically complex production process. Despite state-of-the-art production know-how, the presence of minor spots, air bubbles or streaks in the acrylic panels is unavoidable. In addition, hairline cracks can occur on the acrylic panel during operation.

These occurrences are material-related and are unavoidable in processing. However, they have no significant effect on use and are therefore not accepted as defects.

Over the course of time the acrylic panel will become discolored as a result of usage and must be replaced by customer-service personnel on a regular basis.

The partition works as a UV filter and must be replaced by customer-service personnel on a regular basis.



12.15

How to clean/disinfect

WARNING!



Infections can be transferred through skin contact.

All items/equipment components, with which the user can come into contact during tanning, must be cleaned and disinfected after every tanning session:

- Protective goggles
- Acrylic panels
- Handles and control panel
- AUX cable
- Adjustable air nozzles
- Use only EPA registered quaternary ammonium germicidal based detergent solution (follow product label for use-dilution), such as Lucasol One Step, to clean and disinfect the tanning system and its components.



NOTE:

For thorough cleaning and disinfection, the product must be allowed to remain in contact with the surface being cleaned and disinfected for a specified amount of time.

This specified amount of time depends on the cleaner/disinfectant which is used. Follow the directions for use of the cleaner/disinfectant.

After this specified application time remove remaining cleaner/disinfectant from the surfaces of the device. Use a single-pass paper towel.

Once the remaining cleaner/disinfectant has been removed, perform a visual inspection for signs of any remaining cleaning solution, or areas still soiled. If either are remaining, repeat the cleaning process until a visual inspection can be successfully completed without the need for additional cleaning.

If any of the components listed above cannot be successfully cleaned and disinfected, they must be replaced.



12.15.1 Cleaning and disinfection of surfaces

Acrylic panel

NOTICE

Do not rub with a dry cloth – danger of scratching!

- Other cleaning agents, especially concentrated disinfectants or solvents (e.g. Lysoform, ethyl alcohol or other liquids containing alcohol) must never be used for cleaning. Failure to observe this voids all warranty claims.

Acrylic panels (9) & (11) & (16) referenced in section 12.15.2 'Cleaning/disinfection plan'

These panels are not in direct contact with the user, however they can become soiled if the user inadvertently touches them. If these panels are left soiled for an extended period of time, the heat and light from the tanning lamps may cause the soiled areas to become permanently stained. Therefore we recommend frequently inspecting and cleaning these panels when necessary, to ensure the serviceable life of the acrylic.

- Remove fingerprints and other soiling carefully.

Plastic surfaces

For cleaning the remaining painted and unpainted plastic surfaces, it is best to simply use warm water and commercially available soft chamois leather. Never use aggressive cleaning agents containing alcohol or essential oils. These may lead, in the long run, to damages not covered by the warranty. During cleaning a slight discoloration of the leather cloth by the rubber seals should be expected for production-related reasons.

Prevent damage to the acrylic glass and plastic surfaces.


- Remove rings, watches, bracelets etc. before starting cleaning work.



Filter panels and lamps

- Low pressure lamps can be cleaned with a damp cloth and without additional cleaning agents where necessary.
- High pressure lamps can be cleaned with ethyl alcohol where necessary.

12.15.2 Cleaning/disinfection plan

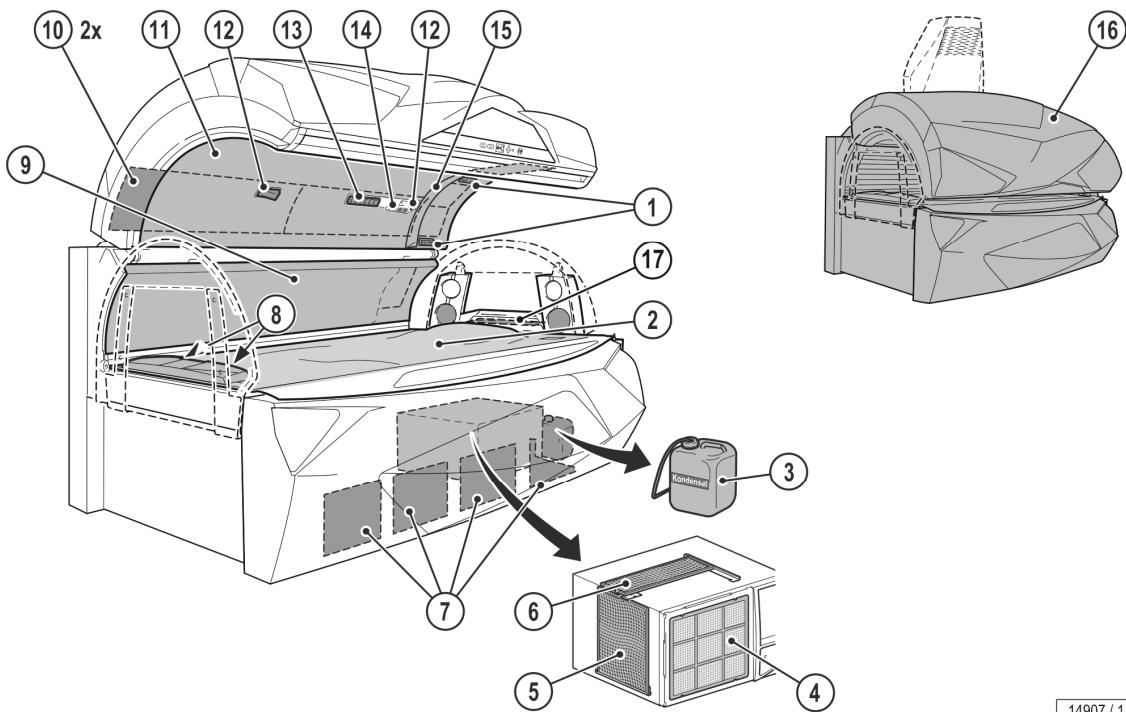

14907 / 1

Component	Interval	Activity	See page
	After each tanning session		
1	X	Clean and disinfect	66
2	X	Clean and disinfect	67
14	X	Clean and disinfect	66
15	X	Clean and disinfect	66
Control Center	X	Clean and disinfect	66
Protective goggles	X	Clean and disinfect	66
AUX cable	X	Clean and disinfect	66
3	X	Empty	90
8	X	Clean and disinfect	66
9, 11, 16	X	Clean	67
Surfaces			
12	X	Clean	66
13	X	Clean	66
4-6	X	Remove and clean filter	90
7	X	Remove and clean filter	88
10	X	Remove and clean filter	89
17	X	Clean	95

**NOTE:**

If desired, you may set the unit to require that cleaning/disinfection be confirmed after each tanning session by pressing the START/STOP button (when the cool-down time for the main fan has finished).



12.16 How to maintain

12.16.1 Technical modifications

The device has been constructed in accordance with current safety regulations. The illustrations and specifications on these instructions are subject to technical modifications which are required in order to improve the device.

Unauthorized modifications and changes to the device and its control unit are prohibited.

Use only OEM spare parts.

12.16.2 Maintenance plan

**NOTE:**

The service of the device may only be carried out by a (Level 2) JK Factory Certified Technician.

NOTICE


Only use OEM spare parts of the same type! The use of other parts will void the manufacturer's warranty!

Any liability is excluded in case of damage or injury proven to be the result of the use of non-genuine spare parts.

The following maintenance activities shall be performed by the studio operator or a (Level 1) JK Factory Certified Technician.

The following items are considered consumables and should be replaced as noted. Replacement of these components helps ensure continued performance of the device.

Maintenance performed by the studio operator or a JK Factory Certified Technician

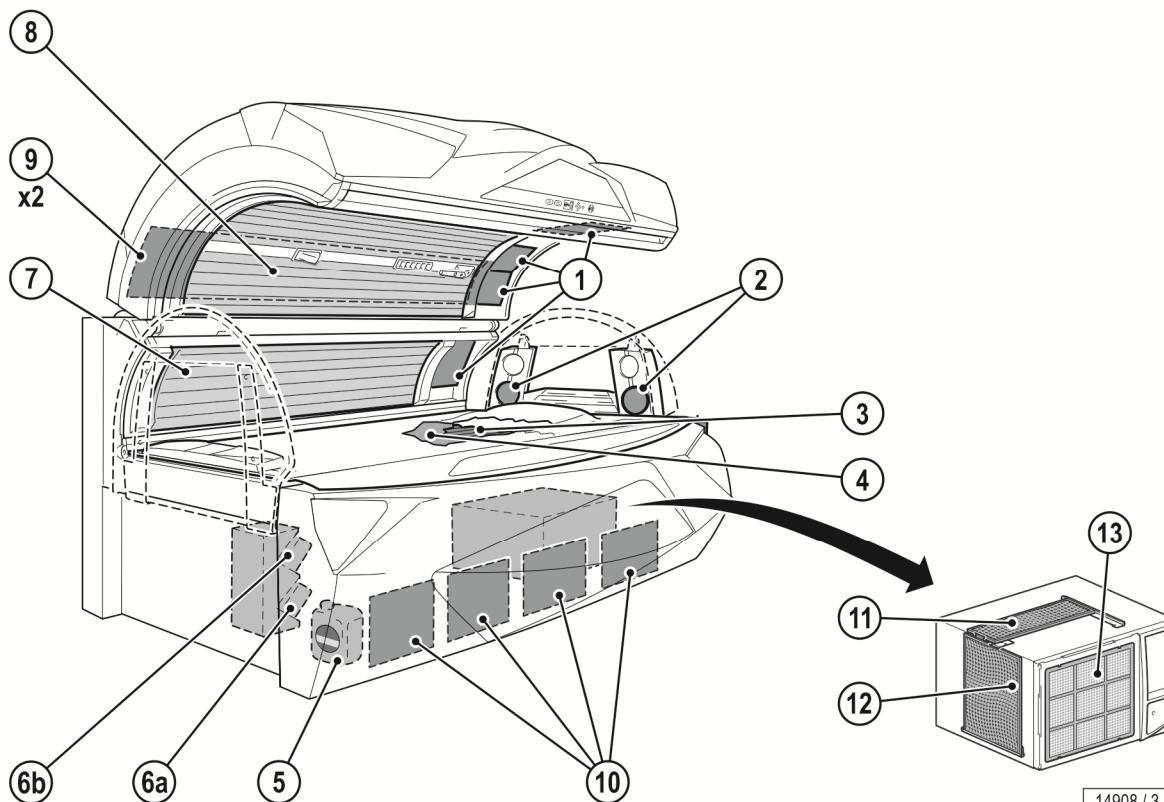

14908 / 3

Component	Interval	When necessary					Activity	See page
		300 Operating hours	500 Operating hours	750 Operating hours	1000 Operating hours	1500 Operating hours		
5 AQUA SYSTEM	X						Replace	93, 94
6a-6b AROMA canister	X						Replace	93
2 Shoulder tanner high pressure lamps		X					Replace	86
1 Facial tanner low pressure lamps		X					Replace	84
1 Facial tanner high pressure lamps			X				Replace	83, 84
3 Low pressure lamps, base				X			Replace	82
7 Low pressure lamps, side section				X			Replace	83
8 Low pressure lamps, canopy				X			Replace	84
1 Starter, facial tanner low pressure lamps					X		Replace	84
1 Facial tanner filter panels						X	Visual control	83, 84
2 Shoulder tanner filter panels						X	Visual control	86

Component	Interval	Activity	See page
	When necessary		
4 Partition	300 Operating hours	X Replace	73
9 Filter mats	500 Operating hours	X Replace	73
10 Filter mats	750 Operating hours	X Replace	73
11 Filter	1000 Operating hours	X Replace	73
12 Filter	1500 Operating hours	X Replace	73
13 Filter	3000 Operating hours	X Replace	73

Maintenance performed by authorized, trained and qualified staff

⚠ WARNING!



WARNING!

Danger to persons from maintenance work that has not been carried out!

If maintenance work is not carried out, it can lead to material damage or personal injuries!

- Maintenance work may only be carried out by personnel who have been authorized, trained and qualified by JK Products & Services, Inc.!
- The operator is responsible for compliance with the maintenance intervals!

The maintenance work must be carried out in accordance with checklist 1019858-.. by a JK Factory Certified Technician at the manufacturer or by personnel who have been authorized, trained and qualified by the manufacturer.

The correctly carried out maintenance work must be entered and confirmed in checklist 1019858-...

The maintenance work must be carried out at the following intervals:

- Every 1500 operating hours
- Annually or every 2000 operating hours, whichever is sooner.
- After 5000 operating hours
- After 10 years or 20,000 operating hours, whichever is sooner.

Once the age of the device reaches ten (10) years from the date of manufacture, its condition must be assessed by a (level 4) JK Factory Certified Technician before use of the device can be continued.


Failure to have this device assessment performed may result in personal injury or property damage.

Any liability is excluded if use of the device continues ten (10) years after the date of manufacture, without the "condition assessment" having been performed.




12.16.3 Cleaning and maintenance work preparation

Service kit

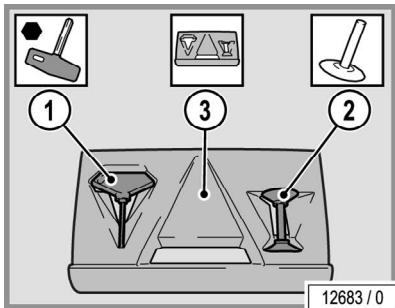


1. Allen key
2. Suction cup for filter panels, acrylic panel and canopy panel
3. Scratch prevention for the entry area (support cushion for the canopy)

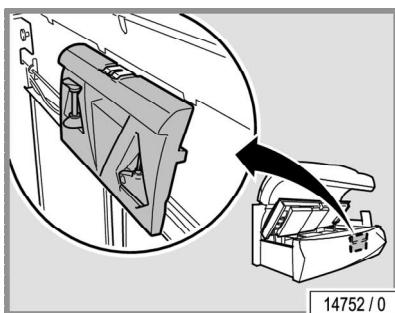


The service kit is fitted on the inside of the front panel.


- Opening the bed base: see page 76.

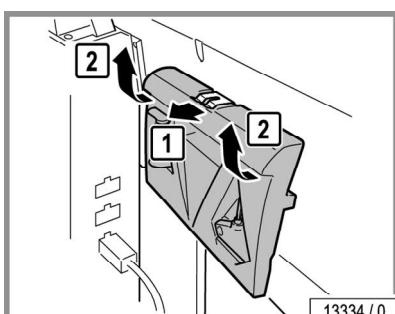


- Remove service kit from support.

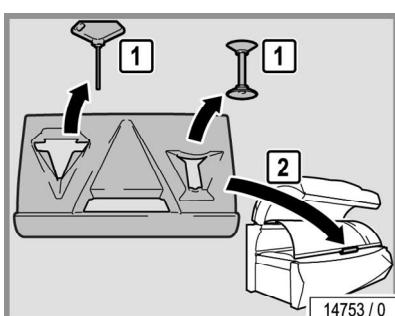

- Remove the Allen-key and suction cup.

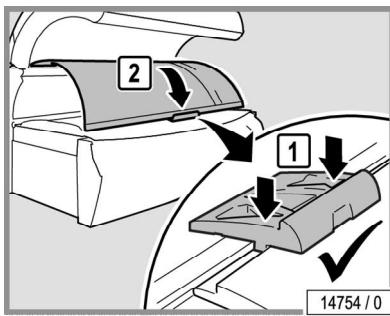


- Use the cushion as scratch prevention when opening the canopy panel.

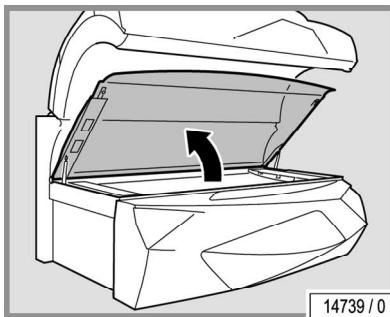


Lift up the acrylic panel and remove the partition


- Lift up the acrylic panel.
The acrylic panel will be held in place automatically by the gas springs.

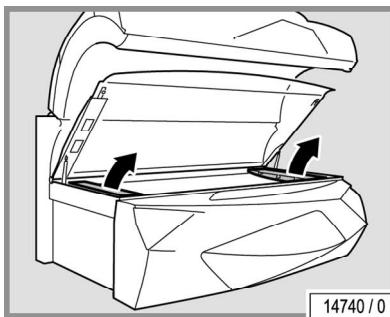

- Remove the lamp covers at the head and foot ends.

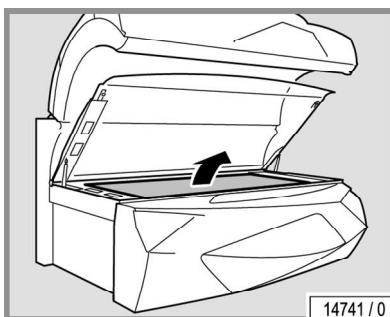


- Remove the partition.

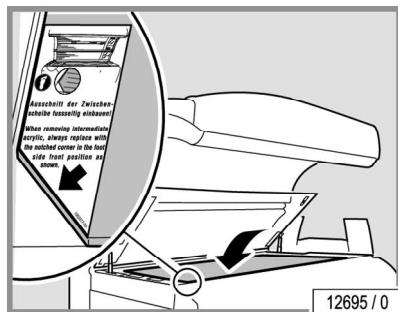


The assembly is the same process in reverse order.

When assembling the device please note:

- You must ensure that the recesses are at the foot end when you fit the partition.

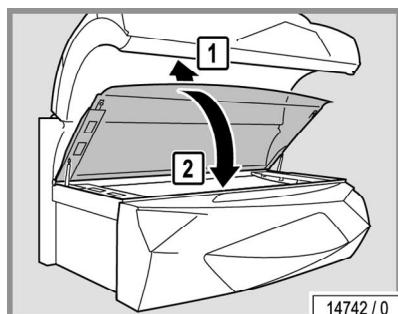


- Lift the acrylic up slightly and then close it.

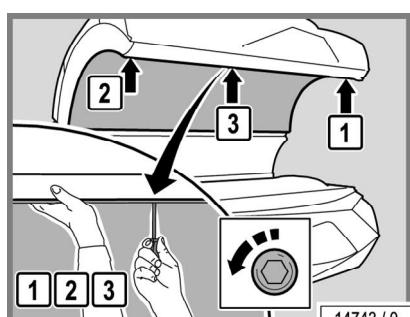


Remove the canopy panel

- Turn the locks through 90°: first the outer locks, then the lock in the middle.
- Support the canopy panel with one hand so that it does not fall down.

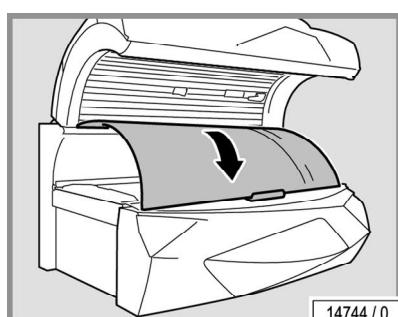

- Carefully lay the panel aside on the scratch protection (service kit).

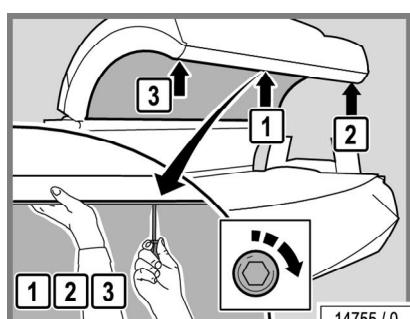

The assembly is the same process in reverse order.



When assembling the device please note:

- Support the canopy panel with one hand so that it does not fall down.
- Re-lock the locks: first the lock in the middle, then the outer locks.




Remove the side section panel

- Turn the locks through 90°.

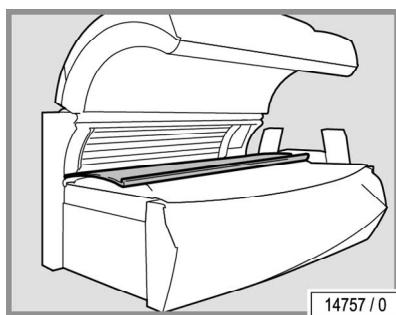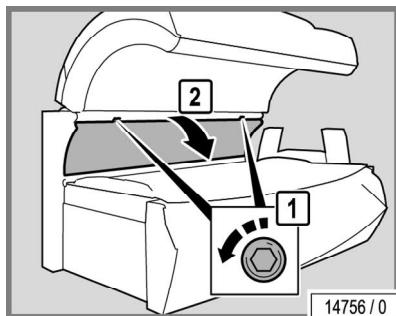

- Open the panel and lay it down carefully.

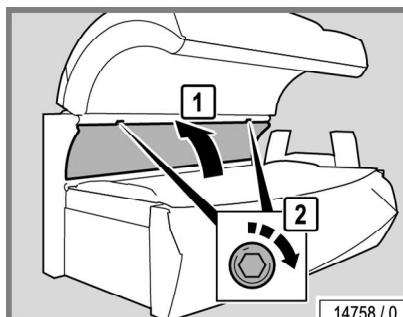


- Re-lock the locks 90°.




Open the bed base


- Press the bed base gently downwards.
- Release the lock at the head end.

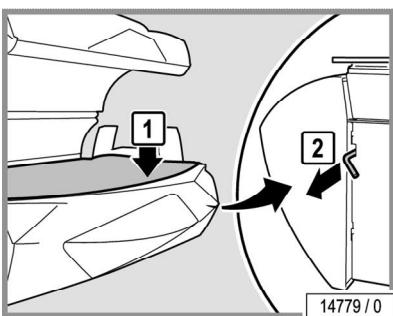

- Lift up the bed base.

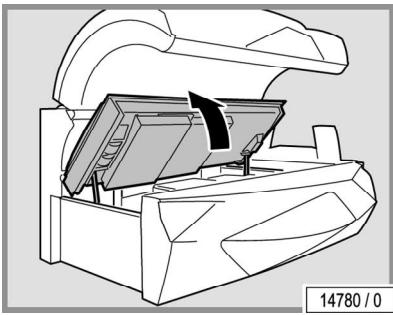

- Close the bed base. You should hear the bed base click into place.

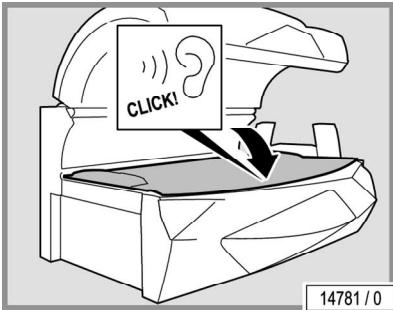



12.16.4 Instructions for replacing the lamps

WARNING!



Danger of being burned by an over-powered lamp!
If the old chip card is left in the device, the new lamp set will be run at too high power and the user may suffer heavy burns!

- The enclosed chip card must always be used when changing the low pressure lamps!

The lamp's UV effect and therefore the tanning effect decrease over time. In electronically controlled low pressure lamps, this loss of effect can be compensated. The chip card is used to control the changing periods. You must adhere to the replacement intervals specified in the maintenance plan.

UV filters are used in order to control the UV radiation, e.g. the filter panels in the facial tanner. However, the filter effect decreases over the course of time. You must adhere to the replacement intervals specified in the maintenance plan.



NOTE:

Small quantities of mercury can escape if the bulb/lamp glass breaks! Please comply with the following instructions:

- Immediately disconnect the device from the mains power.
- Children must leave the room immediately.
- Ventilate the room for at least 15 minutes.
- Put on disposable gloves.
- With the window open: Collect the glass shards and escaped mercury in a sealable container.
- Take the container to a collection site for waste lamps.





Changing the chip cards

The chip card belongs to the long low pressure lamps fitted in the canopy, the side panels and the device's base.

The chip card is located in the 'Part support control'.

- Opening the bed base: see page 76.

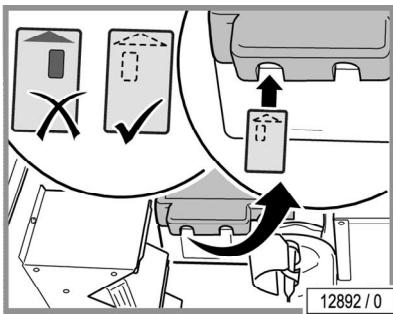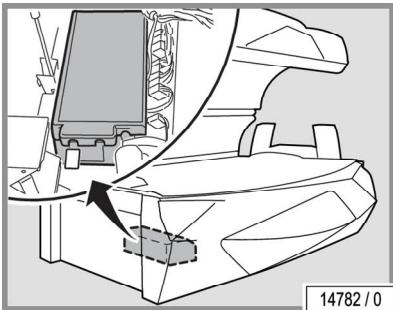

- Remove the chip card belonging to the old lamp set.
- Insert the new chip card with the text facing downwards and the arrow pointing forwards.
- Close the bed base. You should hear the bed base click into place.



Low pressure lamps

⚠ WARNING!







WARNING!

Danger of burns from hot lamps and device parts!

If the tanning device was operated before being switched off and opened, lamps and housing parts may be hot.

- Before replacing lamps or other spare parts, wait until lamps and metal parts of the device have cooled down sufficiently.

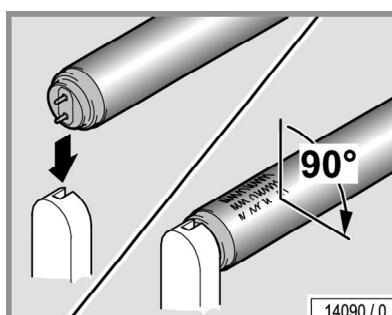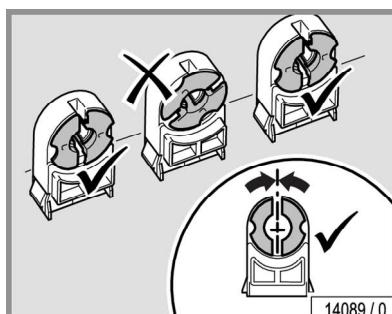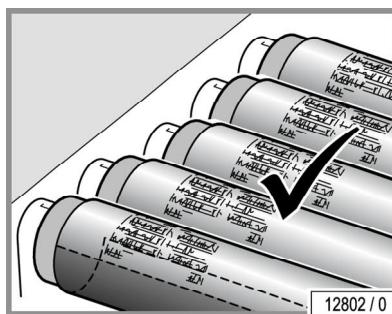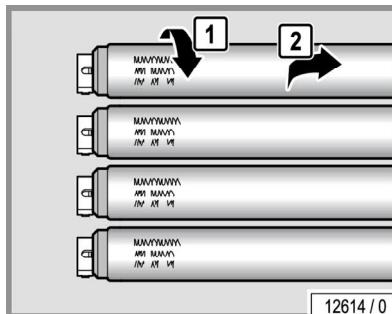

- Replace faulty low pressure lamps immediately.
- Always replace the entire set of lamps after the specified operating time.
- Turn 90° and carefully remove from the socket.

- Low-pressure lamps possess a reflector which directs the light outwards. The lamps must be mounted in such a way that the rear side of the reflector is located on the inside of the unit.

- Push the pins vertically into the socket from above.

- Turn the lamp 90°. The lettering must be oriented towards the outside (i.e. towards the tanning room or towards the user).



Makrolon tubes

Not applicable for this device.

Starter


- Replace faulty starters immediately.

You can recognize a faulty starter by the blackened ends of the lamps or when the lamps can no longer be started.

Operating a lamp with a faulty starter significantly shortens its service life. The lamp loses power and may begin to flicker under certain circumstances.

The starters are located in either one of the two lamp sockets or are visible when the lamps are removed.

- Turn 90° and carefully remove from the socket.

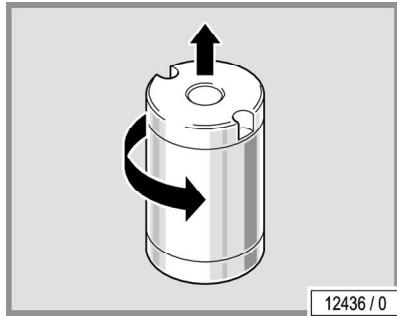



High pressure lamps and filter panels

⚠ WARNING!




WARNING!

Danger of burns from hot lamps and device parts!

If the tanning device was operated before being switched off and opened, lamps and housing parts may be hot.

- Before replacing lamps or other spare parts, wait until lamps and metal parts of the device have cooled down sufficiently.

High pressure lamps are used in the facial tanner and the shoulder tanner.

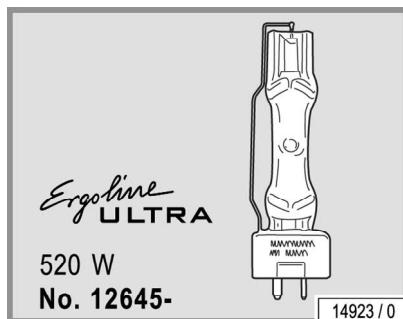

⚠ WARNING!




WARNING!

Fire hazard due to overheated lamps!

Non-approved high pressure lamps can explode. Hot lamp components can set other parts on fire. People can be killed or severely injured due to smoke inhalation and fire.

- Only fit the high pressure lamps specified by the manufacturer.
- Clean the inside of the device regularly. Dust particles are flammable!

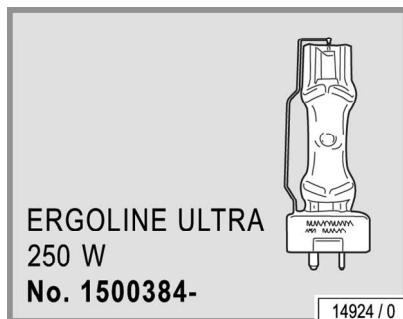


NOTICE

Danger of overheating will exist if the wrong UV high pressure lamps are fitted in the shoulder tanner.


The device might be damaged.

- Never interchange the UV high pressure lamps used in the facial tanner with those used in the shoulder tanner!

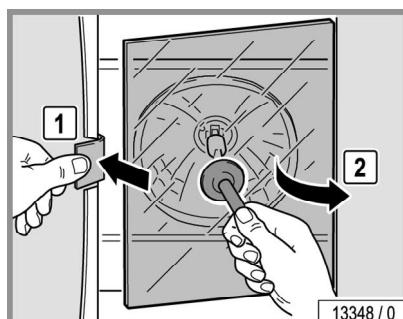


- Release the clamping holder.


- Lift the filter panel, pull it out of the guide and remove it.

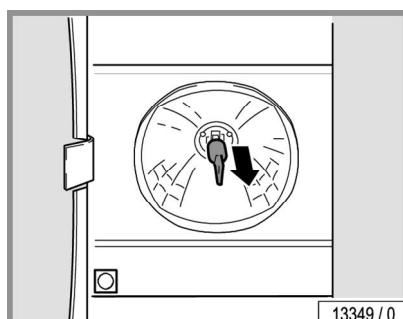


- Unplug the lamps.
- Clean the reflector.

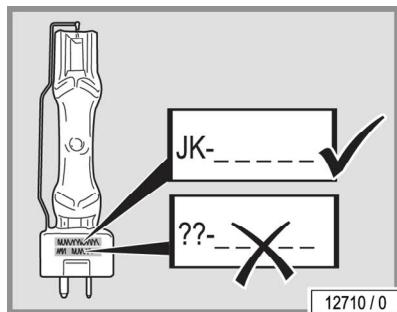

- Only utilize original lamps.

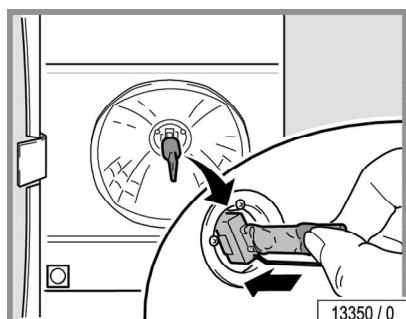

- Only touch new lamps in the handling zone.

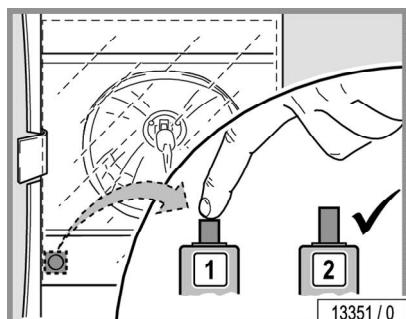

- Test the safety switch to ensure that it works.
The safety switch ensures that the device switches off in the event of a faulty filter panel or when a filter panel is missing.
- Press the button.
If the button springs back automatically then it is functional.
- Replace the filter panel.

12.16.5 Cleaning and disinfection of the sensor


Not applicable for this device.


12.16.6 Cleaning or replacing lamps in the base

- Lift up the acrylic panel and then remove the lamp covers and the partition: see page 73.

Replacing lamps: see page 79.


Changing the chip cards: see page 78.

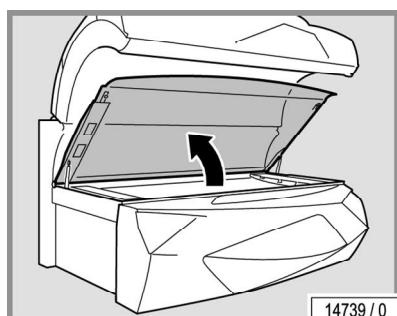

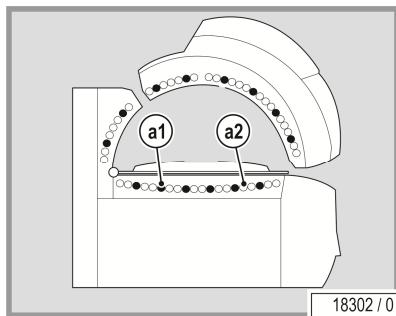

- Please maintain the correct order of the Genesis Type R HYBRID Performance low pressure lamps and Genesis Type U HYBRID Performance low pressure lamps.


a1: Genesis Type R HYBRID Performance low pressure lamps
a2: Genesis Type U HYBRID Performance low pressure lamps

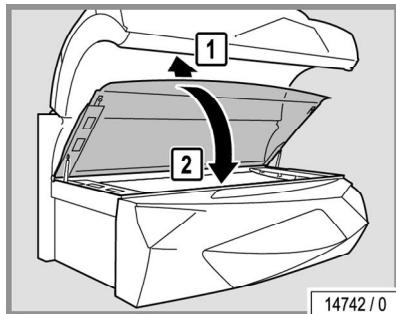

The assembly is the same process in reverse order.


When assembling the device please note:

- You must ensure that the recesses are at the foot end when you fit the partition.
- Lift the acrylic up slightly and then close it.

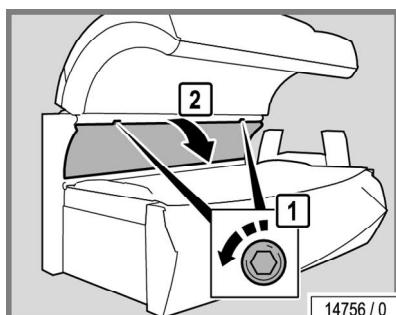

12.16.7 Cleaning or replacing the lamps in the side section

- Removing the side section panel: see page 75.

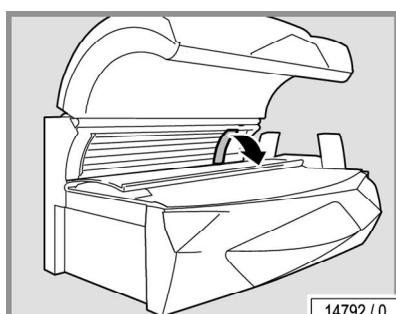

- Remove the lamp cover.


Replacing lamps:

Low pressure lamps: see page 79.


High pressure lamps and filter panels: see page 81.

Replacing lamps: see page 79.

Changing the chip cards: see page 78.

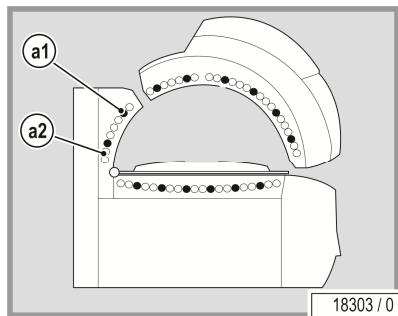


- Please maintain the correct order of the Genesis Type R HYBRID Performance low pressure lamps and Genesis Type U HYBRID Performance low pressure lamps.
a1: Genesis Type R HYBRID Performance low pressure lamps
a2: Genesis Type U HYBRID Performance low pressure lamps

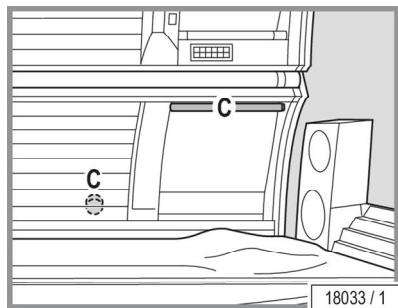

The facial tanner's low pressure lamp must be removed first so that the high pressure lamp can be changed.

- Remove the lamp.

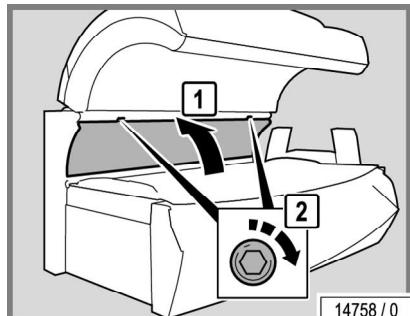

When assembling the device please note:

- Replace and retighten the screws.

12.16.8 Cleaning or replacing lamps and filter panels in the canopy

- Removing the canopy panel: see page 74.

The following components listed in the maintenance plan are fitted in the canopy:


- Low pressure lamps 71 in. (1,80 m)
- Facial tanner low pressure lamps and the associated starters
- Facial tanner high pressure lamps
- Facial tanner filter panels

Replacing lamps:

Low pressure lamps: see page 79.


High pressure lamps and filter panels: see page 81.

Changing the chip cards: see page 78.

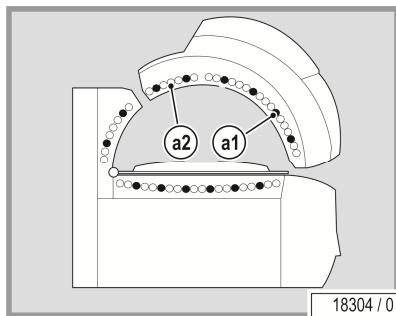


- Please maintain the correct order of the Genesis Type R HYBRID Performance low pressure lamps and Genesis Type U HYBRID Performance low pressure lamps.

a1: Genesis Type R HYBRID Performance low pressure lamps
 a2: Genesis Type U HYBRID Performance low pressure lamps

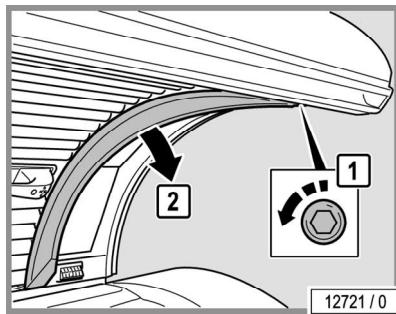


- Remove the screws and then remove the lamp cover.

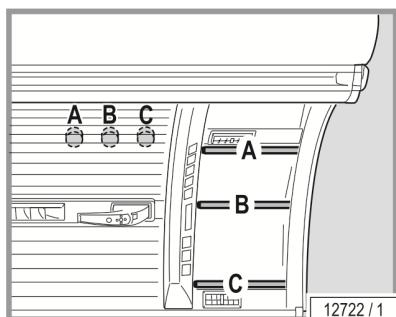


The facial tanner's low pressure lamps must be removed first so that the high pressure lamps and the filter panels can be changed.

- Remove the lamps.

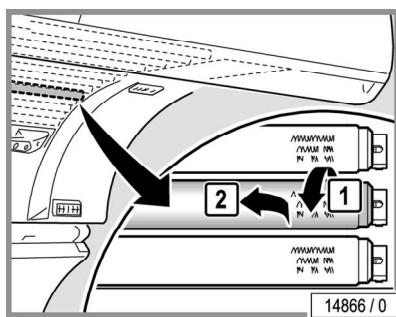

- Remove the ninth low pressure lamp in the canopy so that you can then access the facial tanner's low pressure lamp starter.

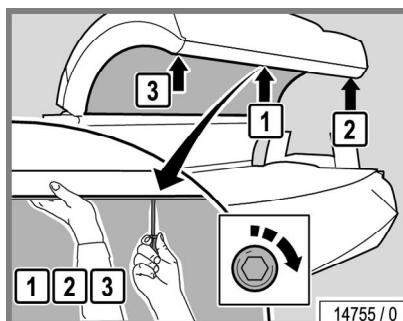

The assembly is the same process in reverse order.

When assembling the device please note:

- Support the canopy panel with one hand so that it does not fall down.
- Re-lock the locks: first the lock in the middle, then the outer locks.

12.16.9 Cleaning or replacing the lamps and filter panels in the shoulder tanner


⚠ WARNING!




Danger of burns from hot lamps and device parts!

If the tanning device was operated before being switched off and opened, lamps and housing parts may be hot.

- Before replacing lamps or other spare parts, wait until lamps and metal parts of the device have cooled down sufficiently.

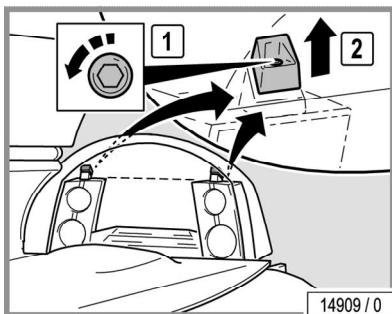

- If an enclosed tanning tunnel is fitted: Remove the screws and then remove the tanning tunnel.

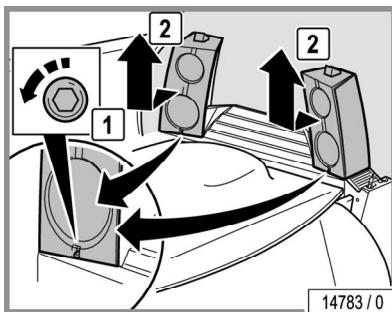

- Remove one screw from each housing.
- Remove both housings and carefully lay them aside.




Changing the high pressure lamps: see page 81.

NOTICE

Danger of overheating will exist if the wrong UV high pressure lamps are used!


The danger of overheating will exist if the wrong UV high pressure lamps are fitted in the shoulder tanner. The device might be damaged.

- Only Ergoline Ultra UV high pressure lamps (1500384-..) are to be fitted in the shoulder tanner.

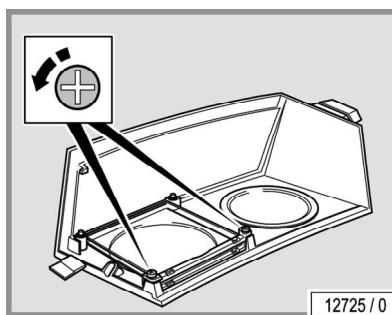


Each of the high pressure lamps in the shoulder tanner is fitted with 2 different filter panels. The two filter panels are secured in place inside the housing.

- Remove the screws.

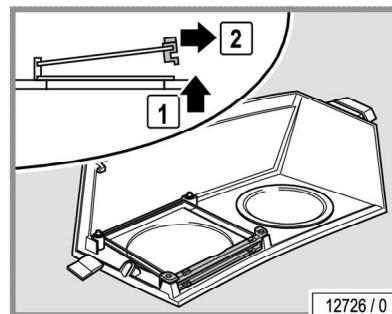

- Use the suction cup to lift the upper filter panel.
- Pull out the holder.

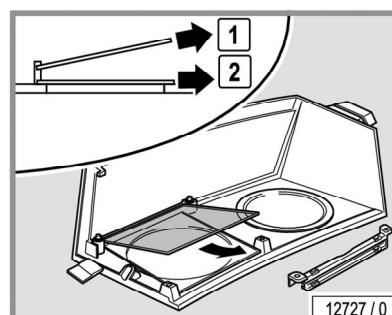


- Use the suction cup to remove the two filter panels.

NOTICE


The shoulder tanner is fitted with different filter panels.

- Ultra Performance 912 (blue), order no. 1500635-..
- Non-coated, order no. 1500160-..

- Do not interchange filter panels!

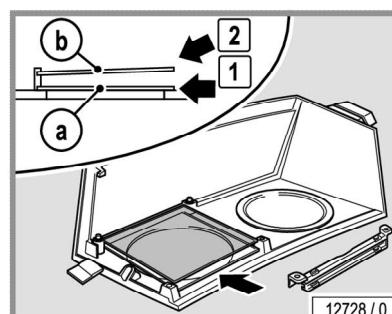

- Fit the new filter panels in the correct sequence.

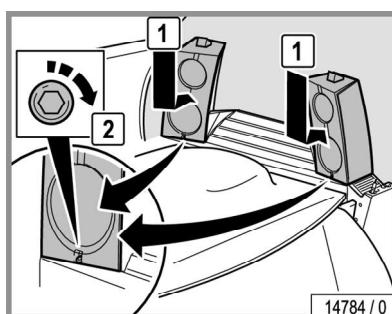

The assembly is the same process in reverse order.



12.16.10 Cleaning or replacing the interior lights

Not applicable for this device.

12.16.11 Cleaning or replacing the canopy accent lighting


Not applicable for this device.

12.16.12 Cleaning or replacing the front panel accent lighting


Not applicable for this device.

12.16.13 Cleaning the filters in the base


- Opening the bed base: see page 76.

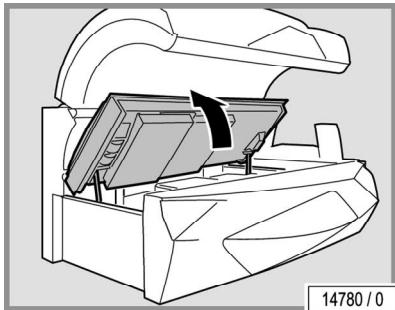


- Remove the filters.

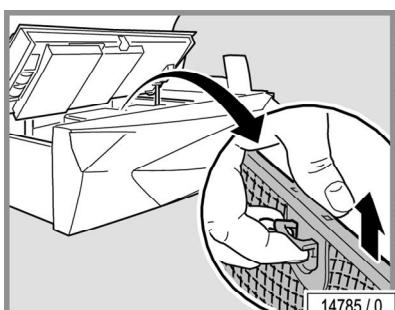

- Clean the filters wet or dry.
Dry cleaning: Vacuum cleaner (depending on the degree of soiling)
Wet cleaning: Use water and detergent or put in a dish-washer

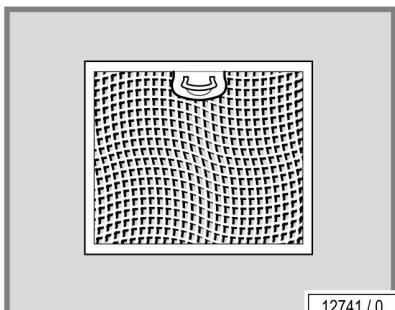

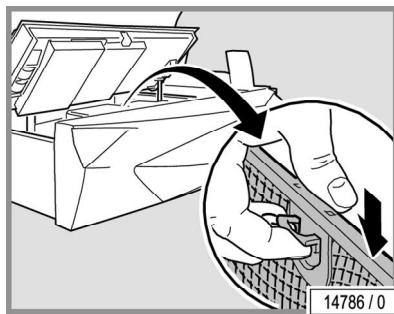


NOTICE!

The device may be damaged by moisture!
The cleaned filters must be dry before they are replaced.




- Reinstall the cleaned and dried filters.
- Close the bed base. You should hear the bed base click into place.

12.16.14 Cleaning the canopy filter mats

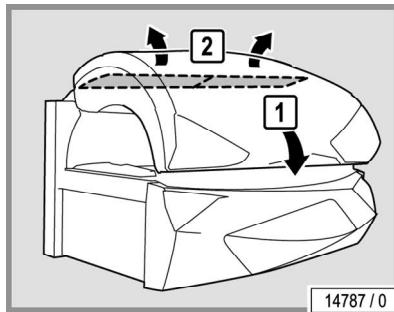


- If necessary, cover the canopy in order to prevent the paint from being scratched.
- Remove the filter mats fitted at the back of the canopy.

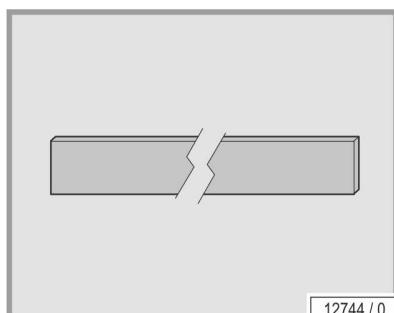

- Dry or moist clean the filter mats or replace with new ones.

Dry cleaning: Vacuum cleaner (depending on the degree of soiling)
Moist cleaning: Water and detergent

- Do not clean in the dish-washer!

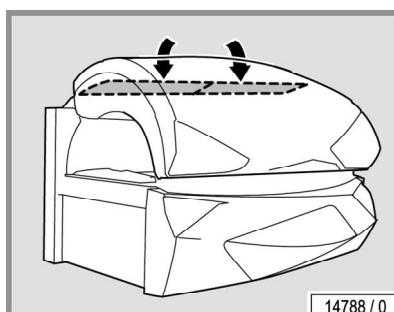

NOTICE!

The device may be damaged by moisture!
The cleaned filters must be dry before they are replaced.

- Refit the cleaned and dried filter mats.


12.16.15 Cleaning the filter mats in shoulder tanner

Not applicable for this device.


12.16.16 Air conditioning: Empty the condensation tank and clean the filter

NOTICE


Damage to the device through soiled filters and cooling fins!
Regularly check the filter grill and cooling fins on the Air conditioning system and clean if necessary.
We recommend using the specially manufactured fin comb, art-no. 15720470, to clean the cooling fins.

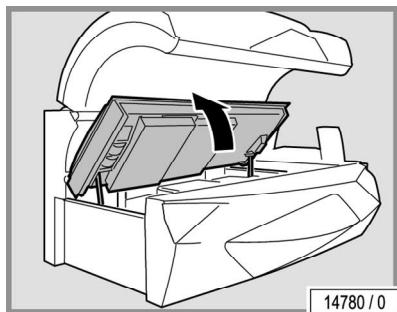


- Opening the bed base: see page 76.

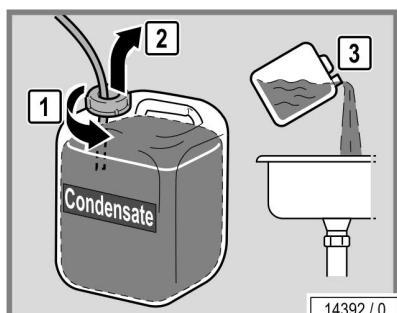


- Bring the condensation tank forwards.
- Disconnect the hose on the condensation tank.
- Empty the condensation tank.

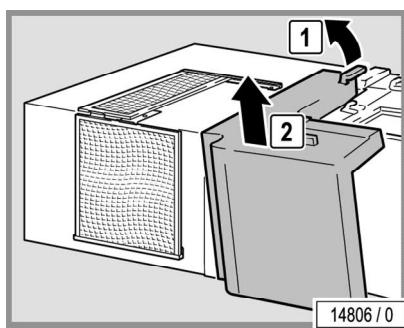


- Remove the cover.

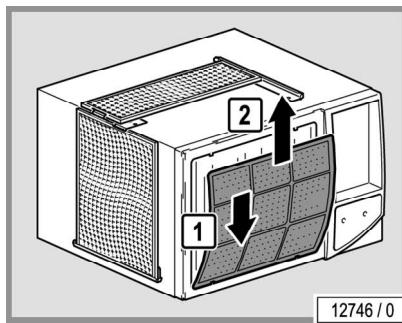


- Remove the filter from the holder.

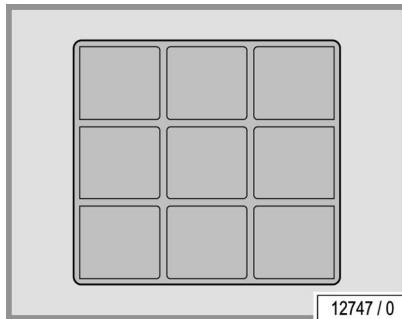


- Clean the filter moist or dry:
Dry cleaning: Vacuum cleaner (depending on the degree of soiling)
Moist cleaning: Water and detergent
- Do not clean in the dish-washer!

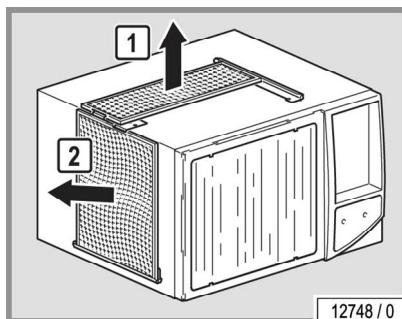

- Remove the filter from the holder.

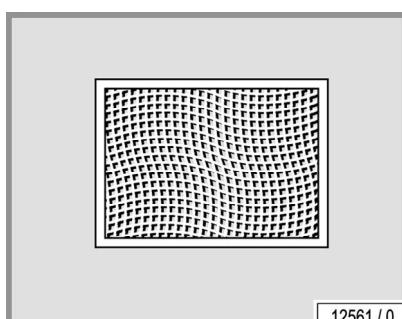


- Clean the filters wet or dry.
Dry cleaning: Vacuum cleaner (depending on the degree of soiling)
Wet cleaning: Use water and detergent or put in a dish-washer

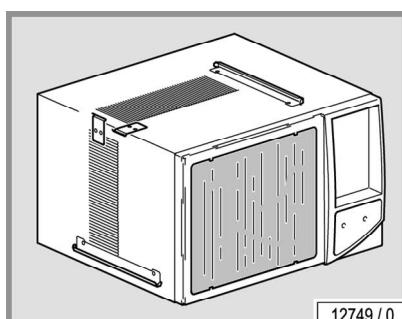

- Use a fin comb to clean the Air conditioning's dirty cooling fins.




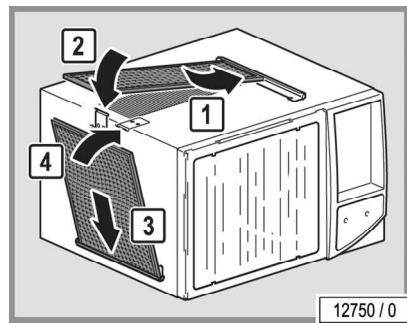


NOTICE!

The device may be damaged by moisture!
The cleaned filters must be dry before they are replaced.


- Reinstall the cleaned and dried filters.
- Ensure that the filters are correctly positioned in the holder.

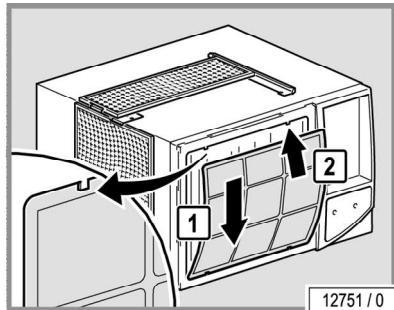

- Attach the cover.

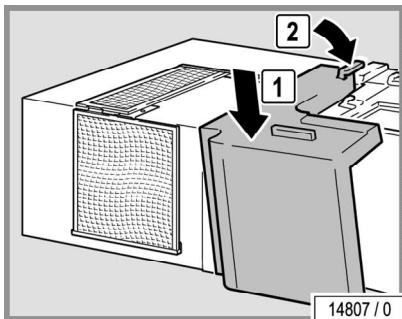

- Connect the condensation tank to the hose.
- Replace the condensation tank.
- Close the bed base. You should hear the bed base click into place.

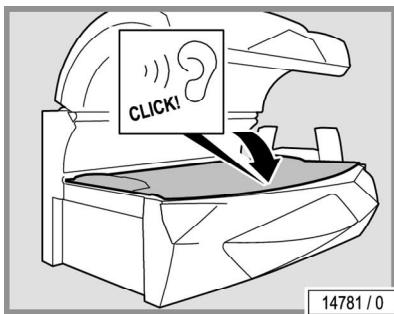



12.16.17 AQUA / AROMA SYSTEM maintenance

The AQUA SYSTEM function sprays a fine mist of AQUA SOLAR WITH SKIN+ ADDITIVE into the interior of the device for body cooling (see 'AQUA SYSTEM PRODUCT INFORMATION', order no. 1011201-..).



⚠ WARNING!



WARNING!

Customer health hazard!



Water that has been contaminated with bacteria can cause skin rashes and other allergic reactions in device users.

- Therefore you must never fill the AQUA SYSTEM canister with tap water or other liquids.
- Only replace the AQUA SYSTEM canister with a genuine canister.
- Do not add any liquid that remains in the old canister to the new canister!
- Ensure that the maximum temperature specified is not exceeded during transport and storage.
- Note the expiration date on the canister.


This AQUA / AROMA SYSTEM also includes the 'RELAX' fragrance, which is distributed through the head nozzles. Another fragrance called 'CABIN' is distributed within the cubicle.

Changing the AROMA canister


- Opening the bed base: see page 76.

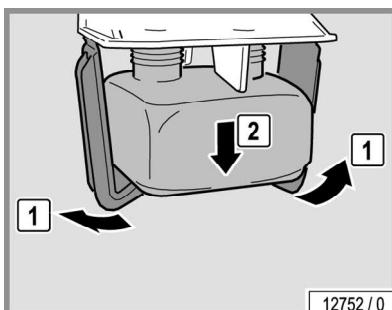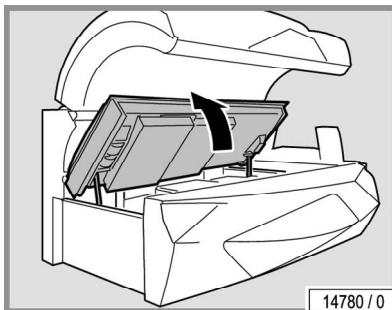


- Remove the AROMA canister.

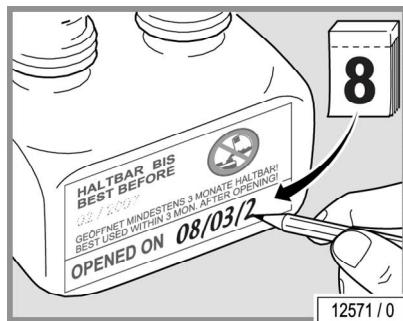

- Unscrew the lid on the new canister.
- Write the date on the canister.

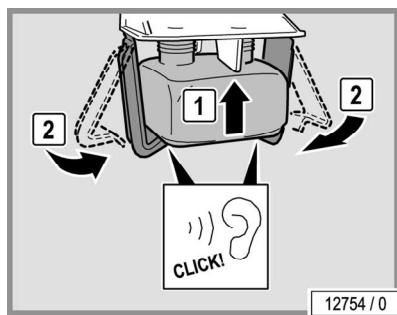


- Fit the new canister in the holder.
- Close the bed base. You should hear the bed base click into place.

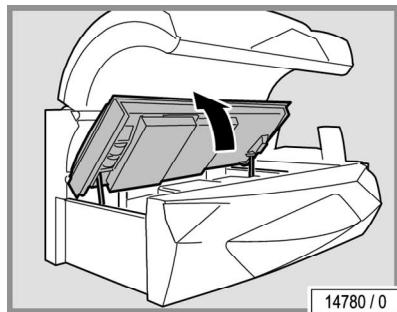


Replacing the AQUA SYSTEM canister

- Opening the bed base: see page 76.

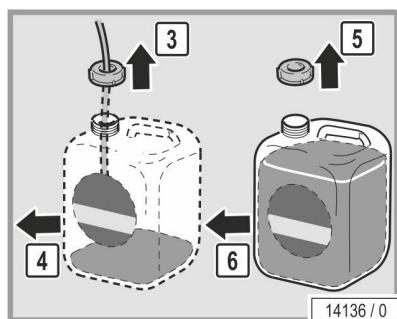

- Remove the canister.
- Unscrew the lid on the new canister.
- Remove the lid with the hose from the old canister.

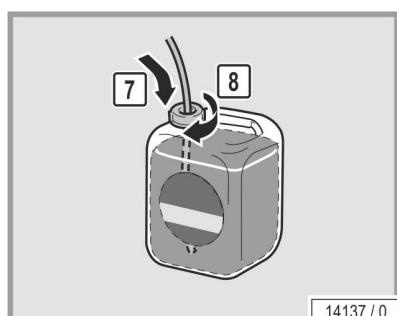

- Screw the lid with the hose onto the new canister.
- Push the hose well down so that it touches the bottom of the canister.
- Put the canister back in the device.
- Close the bed base. You should hear the bed base click into place.

Venting




NOTE:

To ventilate, access to the Professional Setup Manager in the Control Center is required.



- Go to the Professional Setup Manager.
- Activate the fill mode in the AQUA MIST & AROMA menu.

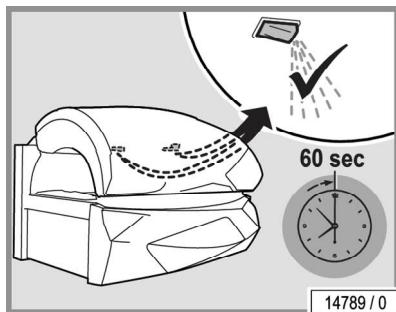


- Close the canopy.

The system will now be ventilated for 60 seconds. The AQUA SYSTEM liquid will be released from the nozzles during this period.

- Then test the AQUA SYSTEM in test mode. In test mode, the nozzles spray for 3 seconds. Repeat the ventilation process if it does not work.

12.16.18 Cleaning the air intake openings

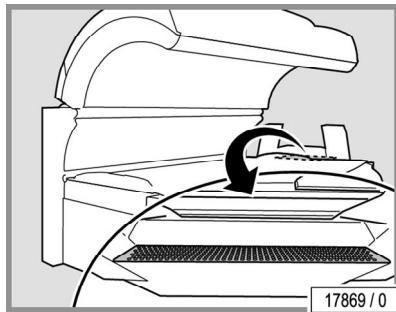


- Use a vacuum cleaner to prevent dust particles.
- Vacuum the air intake openings.



12.17 How to store the device

12.17.1 Decommissioning

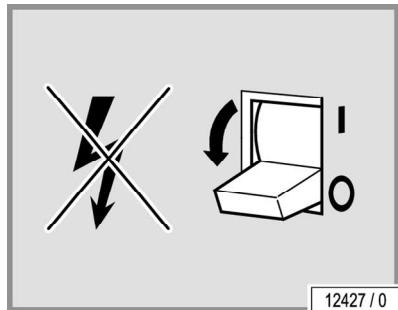

The device must be disconnected from the power supply in order to temporarily or permanently decommission it.

You must abide by the legal disposal requirements when permanently decommissioning the device.

12.17.2 Storage

Store the devices in a dry location with a stable temperature. Once the device has cooled off, pack it in plastic wrapping to protect it against scratches.

- Storage temperature: between 5 °F (-15 °C) and 140 °F (+60 °C)
- Air pressure: between 794 hPa and 1060 hPa
- Relative humidity: between 30 % and 70 %

12.17.3 Disposal

Environmental regulations – disposal of lamps and batteries

Lamps contain fluorescent materials and other waste containing mercury.

Batteries contain heavy metal compounds.

Dispose of lamps and batteries according to the national waste disposal laws and in accordance with the municipal waste regulations.

Packaging

All packaging consists of 100 % recyclable materials.


Disposal of old devices

The device has been manufactured using recyclable materials.

The JK Products & Services will provide you with information on the content or potential hazards of the materials used.

Disposal of consumable supplies

Dispose of consumable supplies according to the national waste disposal laws and in accordance with the municipal waste regulations.

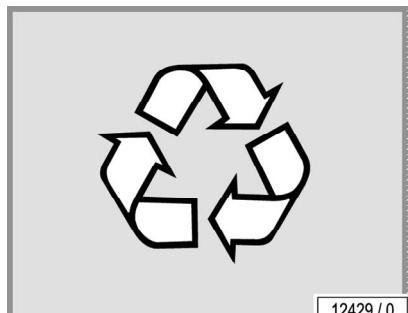



13 What to do when things go wrong

⚠ WARNING!



Follow instructions!

Failure to use protective eyewear may result in severe burns or other eye injury. If discomfort develops, discontinue use and consult a physician.

- A person who does not tan in the sun most likely will not tan from the use of this device. Persons with very sensitive skins (skin type I) must not use the unit, as with this skin type tanning is impossible without damaging the skin (e.g. sunburn).
- Individuals with acute sunburn, individuals who have had or have skin cancer, or individuals who are at increased risk of skin cancer, must not use the system.
- UV devices must not be used without taking medical advice if unexpected effects, such as itching, occur within 48 hours of the first session.
- The tanning device must never be used if the timer is defective!
- UV light from the sun or UV systems can cause injuries to the skin or eyes. This biological effect is dependent on the sensitivity of the individual's skin and on the type and amount of UV light.
- As with natural sunlight, overexposure can cause eye and skin injury, sunburn, and allergic reactions.
- Repeated exposure may cause premature aging of the skin and skin cancer.
- Medications (taken internally or otherwise) or cosmetics applied to the skin, may increase sensitivity to ultraviolet light.
- A person having a history of skin problems should consult a physician before use.

- Seek medical advice if in doubt!



NOTE:

Occasionally, persons using the tanning device will experience a slight reddening of the skin – usually in small patches – after the second or third session. This redness is often accompanied by an itching sensation. This may be nothing more than a very limited heat 'rash' caused by heat from the lamps within the device. It should go away within approximately 24 hours and should not reappear. This rash can be lessened or prevented by applying moisturizer lotion to the affected area after the tanning session is completed.

- Follow instructions.
- Increasing tanning also requires an extension of the exposure time (= radiation time) or after a certain degree of tanning has been achieved, a further intensification is not possible. However, the radiation time may not be extended as desired within the scope of the permissible radiation doses! Therefore, only a certain degree of final tanning dependent on the skin type can be achieved without endangering your health.
- Inform users about proper tanning.



13.1 Troubleshooting

Software and device malfunctions are indicated by error codes, e.g. #0701.

These error codes are output to the display to simplify troubleshooting:

- If a fault occurs, the error code flashes on the display.
- If several errors occur, the error messages are displayed alternately.
- Error elimination is acknowledged by pressing the START/STOP button.
- If the error cannot be cleared, contact Customer Service (see page 2) or your local service provider

More detailed information regarding the error description and solutions is available in the document 'Error codes' #1011831-.., which you find in the Technical Documentation folder.

For authorized, trained and qualified staff:

For further technical description, please contact Customer Service (see page 2).



NOTE:

As a result of EMC phenomena (high-frequency irradiation and surge) the audio system may malfunction.

This can be rectified by switching off the power supply and switching it back on again (power on reset).



WARNING!

Danger of skin injuries and skin diseases.



WARNING!

- The tanning device must never be used if a filter panel is missing or damaged, or the timer is defective!

14 Where you can find out more



If you have any questions, please contact Customer Service (see page 2).

15 How we studied this device



The system is built to the state of the art and the applicable safety regulations.

We reserve the right to make technical changes which are necessary to improve the system or a system component even if they deviate from pictures and information in these instructions.

16 Appendix



16.1 Model-No.: JK 162 / 52-4 TT HP: EMC information

16.1.1 Table 201: Guidelines and manufacturer's declaration – Electromagnetic emissions

The tanning system has been designed for use in the electromagnetic environment specified below. The device operator must ensure that it is used in such an environment.

Emission measurements		Conformity Electromagnetic environment – Guidelines
HF emissions as per CISPR 11	Group 1	Model-No.: JK 162 / 52-4 TT HP uses HF power mostly for internal functions and Bluetooth®. Therefore the HF emissions are very low and it is highly unlikely that it will cause any interference to adjacent electronic equipment.
HF emissions as per CISPR 11	Class B	Model-No.: JK 162 / 52-4 TT HP has been designed for use in all devices, including those in residential areas, and it has been determined that it can be used for residential purposes as soon as it is connected up to a public mains supply that also supplies buildings.



16.1.2

Table 202: Guidelines and manufacturer's declaration – Electromagnetic interference

The tanning system has been designed for use in the electromagnetic environment specified below. The device operator must ensure that it is used in such an environment.

Interference immunity tests	IEC 60601 - Test level	Defined level	Electromagnetic environment – Guidelines
Discharging static electricity (ESD) IEC 61000-4-2	±8 kV contact discharge ±2 kV air discharge ±4 kV air discharge ±8 kV air discharge ±15 kV air discharge	±8 kV contact discharge ±2 kV air discharge ±4 kV air discharge ±8 kV air discharge ±15 kV air discharge	–
Fast transient electrical noise / bursts IEC 61000-4-4	±2 kV for mains supplies	±2 kV for mains supplies	–
Surges 61000-4-5	±1 kV series mode voltage ±2 kV common mode voltage	±1 kV series mode voltage ±2 kV common mode voltage	–
Magnetic field radiated by supply frequency (50/60 Hz) IEC 61000-4-8	30 A/m	30 A/m	–
Voltage drops, Short interruptions and supply voltage IEC 61000-4-11	<5% U_T (>95% drop in U_T) for 5 s Only this test is performed, because the device used 2 or 3 phases and has more than 16A / phase.	<5% U_T (>95% drop in U_T) for 5 s	–

Note: U_T is the mains alternating voltage used for the test levels.



16.1.3 Table 204: Guidelines and manufacturer's declaration – Electromagnetic interference

The tanning system has been designed for use in the electromagnetic environment specified below. The device operator must ensure that it is used in such an environment.

Interference immunity tests	IEC 60601 - Test level	Defined level	Electromagnetic environment – Guidelines
Conducted HF noise IEC 61000-4-6	3 V rms 150 kHz to 80 MHz outside the ISM band, 6 V RMS in the ISM and amateur radio bands	3 V rms 150 kHz to 80 MHz outside the ISM band, 6 V RMS in the ISM and amateur radio bands	Portable and mobile radios are not to be used at shorter distances to the Model-No.: JK 162 / 52-4 TT HP and its cables other than the recommended safe distance, which has been calculated according to the transmitting frequency are made.
	Modulation: 80% AM at 1kHz	Modulation: 80% AM at 1kHz	Recommended safe distance: $d = 1.2\sqrt{P}$
Radiated HF noise IEC 61000-4-3	10 V/m 80 MHz to 2,7GHz	10 V/m 80 MHz to 2,7GHz	$d = 1.2\sqrt{P}$ 80 MHz to 800 MHz
	Modulation: 80% AM at 1kHz	Modulation: 80% AM at 1kHz	$d = 2.3\sqrt{P}$ 800 MHz to 2,7 GHz
	27 V/m 385 MHz	27 V/m 385 MHz	P is the transmitter's maximum output power (nominal rating) in Watts (W) according to the details provided by the manufacturer of the transmitter and the recommended safe distance is given in metres (m).
	Modulation:	Modulation:	The field strength of the stationary radio transmitter is for all frequencies ^b in accordance with tests that were carried out on-site ^a less than the determined level.
	Pulse Modulation: 18 Hz	Pulse Modulation: 18 Hz	Interference is possible in the immediate surrounding of equipment that carry the following symbols.
	28 V/m 450 MHz	28 V/m 450 MHz	
	Modulation:	Modulation:	
	FM + 5 Hz deviation: 1 kHz sine	FM + 5 Hz deviation: 1 kHz sine	
	9 V/m 710, 745, 780 MHz	9 V/m 710, 745, 780 MHz	
	Modulation:	Modulation:	
	Pulse Modulation: 217 Hz	Pulse Modulation: 217 Hz	
	28 V/m 810, 870, 930 MHz	28 V/m 810, 870, 930 MHz	
	Modulation:	Modulation:	
	Pulse Modulation: 18 Hz	Pulse Modulation: 18 Hz	
	28 V/m 1720, 1845 MHz	28 V/m 1720, 1845 MHz	
	Modulation:	Modulation:	
	Pulse Modulation: 217 Hz	Pulse Modulation: 217 Hz	
	28 V/m 1970, 2450 MHz	28 V/m 1970, 2450 MHz	
	Modulation:	Modulation:	
	Pulse Modulation: 217 Hz	Pulse Modulation: 217 Hz	
	9 V/m 5240, 5500, 5785 MHz	9 V/m 5240, 5500, 5785 MHz	
	Modulation:	Modulation:	
	Pulse Modulation: 217 Hz	Pulse Modulation: 217 Hz	

Comment 1: The higher value applies to 80 MHz und 800 MHz.

Comment 2: These guidelines might not apply to all situations. The propagation of electromagnetic waves is affected by absorption and reflexes caused by buildings, objects and people.

a The field strength of a stationary transmitter e.g. standard cell phone stations (mobile / cordless phones) and public mobile radio services, amateur stations, MW and UHF radio and TV transmitters) cannot be precisely determined for theoretical reasons. It is recommended that an onsite test should be carried out in order to determine the electromagnetic environment caused by stationary HF transmitters. If the field strength of the massage device determined at the installation site exceeds the defined HF level listed above then the device must be monitored with regard to its normal operation. If any abnormal characteristics are recorded then it will be necessary to introduce additional measures, such as reorientation or converting the device.

b The field strength must be less than 3 V/m over the frequency ranges from 150 kHz to 80 MHz.



16.1.4

Table 206: Recommended safe distances between portable and mobile HF communications equipment and Model-No.: JK 162 / 52-4 TT HP

Model-No.: JK 162 / 52-4 TT HP has been designed for use in an electromagnetic environment in which the radiated HF interference is controlled. The operator of the device can help to reduce electromagnetic interference by maintaining a safe distance between the portable and mobile HF communication devices (transmitters) and the device itself as recommended below with regard to the communication device's maximum output power.

Nominal transmitter power (W)	Safe distance with regard to the transmitting frequency (m)		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800 MHz to 2.7 GHz $d = 2.3\sqrt{P}$
0,01	0,12	0,12	0,23
10	3,8	3,8	7,3
1	1,2	1,2	2,3
0,1	0,38	0,38	0,73
100	12	12	23

The recommended safe distance (d) in metres (m) can be determined by using the equation that applies to the transmitters frequency in the case of transmitters where their nominal power is not listed above, whereby 'P' is the nominal transmitter power in Watts (W) according to the details supplied by the manufacturer of the transmitter.

Comment 1: The higher value applies to 80 MHz und 800 MHz.

Comment 2: The propagation of electromagnetic waves is affected by absorption and reflexes caused by buildings, objects and people.

EM-Emission inside device

5 Hz – 2 kHz	max. 2.1 µT (1.5 kV/m)
2 kHz – 400 kHz	max. 3.3 µT (1.3 kV/m)
9 kHz – 30 MHz	max. 35 dBµV/m (56 µV/m)
30 MHz – 1 GHz	max. 95 dBµV/m (56 µV/m)
1 GHz – 6 GHz	max. 90 dBµV/m (32 µV/m)

See also page 38.

Influence of electromagnetic disturbances

If the essential performance features (see 3.10.1) are influenced by electromagnetic interferences or do not exist at all, there is no risk for the user.

17 JK-Products, Inc. – Limited Warranty



JK-Products, Inc. warrants its products to be free from defects in materials and workmanship under intended normal use as described in the unit's Operation and Installation Instruction, for a period of three (3) years from date of sale.

This Limited Warranty applies only to the original purchaser of the equipment through JK-Products, Inc. or its authorized dealer or distributor, and is not transferable.

Global Service USA obligations under this warranty are limited to repair or replacement of any defective part without charge for that part to the original purchaser, with the following exceptions:

A Fluorescent lamps are warranted against defects for a period of thirty (30) days from date of sale.

B Only parts obtained through JK-Products, Inc., its authorized dealers or distributors may be used. Transportation costs for parts shipped to the consumer and the return of defective parts to Global Service USA are not included.

C Labor will be furnished without charge for ninety (90) days from date of purchase only. All labor and related charges must be authorized by Global Service USA prior to start of repairs, and must coincide with Global Service USA established rates and time allotment policy.

D Acrylic: Refer to Manufacturer's Acrylic Warranty Policy. It is imperative that the original customer completes and returns the enclosed warranty card within 10 days after purchase to insure valid registration and coverage for potential claims.

E Items considered to be wear items such as rubber foot mats, buttons, or other items which the user comes into contact with during normal operation of the device, are warranted for a period of 30 days from the date of purchase.

If the warranty card is not registered, proof of purchase from JK-Products, Inc. or its authorized dealer or distributor will be required prior to any consideration on warranty claims. This could result in service delays.

This warranty is extended to the individual or legal entity, whose name appears on the warranty registration card filed with Global Service USA, or whose name appears on the original sale document and may not be transferred to any other individual or legal entity.

This warranty does not apply to any failure of the product, or any parts of the product due to alterations, modifications, misuse, abuse, accident, improper maintenance, improper installation, acts of God or if the serial number on the product has been removed, altered or defaced. Adequate packaging must be used for returned goods to prevent freight damage.

There are no warranties which extend beyond the description on the face hereof. There is no express or implied warranty of fitness of merchantability.

The remedies provided in this limited warranty are the exclusive remedies provided to the purchaser by Global Service USA and are provided in substitution of all other remedies. Consequential and incidental damages are excluded.

No person, firm or corporation is authorized to obligate JK-Products, Inc. for any liability in connection with the sale or use of these goods.

17.1 Limited Acrylic Warranty



JK Products, Inc. warrants its acrylic sheets to be free from defects in material and workmanship, under intended normal use, for a period of one (1) year from date of sale of the tanning bed or 1,200 hours of operation whichever comes first. Due to the tanning lotions, cosmetics, disinfectant and improper cleaners used on tanning surfaces that cannot be controlled by JK Products, Inc., after the first thirty (30) days of warranty period acrylic will be subject to a prorated cost to the consumer.

JK Products, Inc. obligations do not include transportation charges on replacement acrylic.

JK Products, Inc. assumes no liability for the cost of removing defective sheets or installing replacement sheets, or for damage to persons or property. To make a claim, contact an authorized dealer or distributor of JK-Products, Inc. or Global Service USA with the model number, the serial number, and the date of purchase of the bed in which the acrylic will be used.

There are no warranties which extend beyond the description on the face hereof. There is no express or implied warranty of fitness or merchantability.

The remedies provided in this limited warranty are the exclusive remedies provided to the purchaser by Global Service USA and are provided in substitution of all other remedies. Consequential and incidental damages are excluded.

JK Products & Services, Inc.
411 West Washington Avenue
Suite B
Jonesboro, Arkansas 72401
Phone: 800-445-0624
Fax: 870-935-9046
E-mail: support@jkamerica.com
<http://jkproducts.us/service>



17.2 Warranty Claims

17.2.1 Policy & Procedures

1. When to file a Warranty Claim:

In the event that your tanning system is not functioning properly under intended normal use as described in the unit's Operation and Instruction Manual, you may then have established cause to file a Warranty Claim.



NOTE:

Both the tanning system and legal registered owner must meet the criteria established under JK Products, Inc. Limited Warranty as described in this manual before any Warranty Claim will be considered.

2. Where to call to place a Warranty Claim:

JK-Products, Inc. is supported nationwide by an extensive dealer/distributor network through whom the majority of our units are sold. These dealers/distributors are required by contract to stock sufficient parts to meet the various warranty and nonwarranty needs of their customers.

Therefore, in order for you to properly file a Warranty Claim, you must locate the name and phone number of the dealer/distributor through whom you purchased your unit. Normally this information can be easily located either on your paid invoice or other proof of purchase documentation that you received when taking delivery of the unit.

Next, notify the dealer/distributor of the serial number, model number and date unit was purchased. This will initiate the Warranty Claim process.

Upon determination of the problem, your dealer/distributor can then supply you with the part(s) and proper instructions to return your unit to working order.

Each dealer/distributor has their own internal procedure for handling Warranty Claims and credits.



NOTE:

JK-Products, Inc., the manufacturer, does not involve itself in the direct handling of a Warranty Claim except where mitigating circumstances apply and where required by law.



NOTE:

If a part returned to Global Service USA proves to be operational and not defective, or if product code numbers on part do not correspond with individual unit date on file with Global Service USA, then a 30% restocking fee will be charged in addition to the retail purchase price of part.

Acknowledgment of the Agreement is not required for acceptance of the foregoing stipulations.

You must comply to the instructions and warnings under chapter:

- 4.1 Intended use
- 5.1 Prohibited usage
- 6. What you must do to avoid serious harm



Any other use shall be considered improper. The manufacturer cannot be held liable for damage or injuries resulting from misuse. The operator bears the sole risk for this.

The proper use also includes compliance with the manufacturer's instructions, operating and maintenance conditions. This device may only be operated, maintained and repaired by persons familiar with these tasks and who are aware of the dangers involved.

17.2.2 Labor Warranty Claim/Credit Procedures


Global Service USA furnishes labor without charge for a period of 90 days from date of purchase. All labor and related charges must be authorized by Global Service USA prior to start of repairs, and must coincide with Global Service USA established rates and time allotment policy.

Please contact our Service Department for authorization and billing procedures.

17.2.3 Quality inspection and acceptance form

The form can be found in the separate document 'Quality inspection and acceptance form' (#1014113-..) which is provided in the Technical Documentation folder.

Tear here to remove page

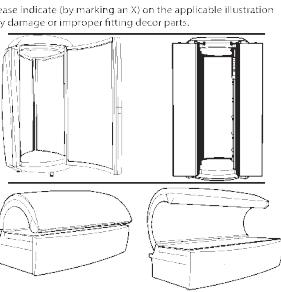Quality Inspection and Acceptance Form. Both parties must sign the form in order for it to be valid. Upon completion, email the signed form to JK Products & Services at Questions? (800) 643-0086 Q.F@JK-AMER-CA.COM		Please indicate here if device rebranded and labels applied: <input type="checkbox"/>
Client/Salon Name: _____ Contact: _____ Street Address: _____ City: _____ State: _____ Zip: _____ Phone Number: _____ Install Date: _____ Service Provider: _____		
Select the box next to the appropriate subsystem indicating the system is installed and functioning properly. If the device is not equipped with a particular subsystem draw a line through it. Lamps: <input type="checkbox"/> Body Bench: <input type="checkbox"/> Body Canopy: <input type="checkbox"/> Facial (HP): <input type="checkbox"/> Facial (Spa/Spa) <input type="checkbox"/> Shoulder Tanner: <input type="checkbox"/> Decor: <input type="checkbox"/> Options: <input type="checkbox"/> Skin Sensor: <input type="checkbox"/> Air Conditioning: <input type="checkbox"/> Aqua Fresh: <input type="checkbox"/> Aroma: <input type="checkbox"/> Audio System Select the box next to the appropriate subsystem indicating the system has been setup and adjusted. If the device is not equipped with the subsystem draw a line through it. <input type="checkbox"/> Lifting Mechanism: <input type="checkbox"/> Timer Network: _____		Affix Name Plate from Manual: HERE:
Select the box next to the operating feature and maintenance functions that have been explained to the customer or salon representative. <input type="checkbox"/> Basic Operation: <input type="checkbox"/> Accessing lamps: <input type="checkbox"/> Air Filter Location: <input type="checkbox"/> Feature Operation: _____ <input type="checkbox"/> The Operating Instructions, Service Tool, and Eye Protection have been provided. <input type="checkbox"/> The device has been installed, electrically connected and fully operational. <input type="checkbox"/> The device has been installed, but not electrically connected all systems appear to be operational.		Acceptance testing: Equipment Used: _____ Serial Number: _____ Calibration Due: _____ High Voltage (Hypot): _____ Ground Bond: _____ Ammeter: _____ Voltmeter: _____ High Voltage (Hypot) Test: 1.25 KV AC / 2.5 KV DC <input type="checkbox"/> Pass <input type="checkbox"/> Fail Ground Bond: Less than 0.1 Ohms <input type="checkbox"/> Pass <input type="checkbox"/> Fail Amperage: Line 1: _____ Amperage: Line 2: _____ Amperage: Line 3: _____ Voltage: L1 to L2: _____ Voltage: L1 to L3: _____ Voltage: L2 to L3: _____
Signing this form, confirms the device was installed and tested both electrically and functionally by the Service Provider, any damage is indicated in the section to the right, and the collections above have been completed, explained, and understood by all. Service Provider Representative Print Name: _____ Signature: _____ Date: _____ Salon/Club Representative Print Name: _____ Signature: _____ Date: _____		

03/2019 1014113-02- Quality Inspection and Acceptance Report

19320 / 0



17.3 Limited Lifetime Warranty

Global Service USA warrants specific individual components of the tanning system against defects in materials and workmanship for the life of the product. Coverage is limited to units purchased after November 1, 1994.

This warranty is restricted to the following components only:

- Steel Back Panels
- Structural Steel Frame Components
- Exterior Steel Body.



NOTE:

Components are warranted on an individual basis and not as they may effect the complete product.

Warranty coverage does not include cosmetic abnormalities such as scratches, nicks, dents, rust, oxidation of paint or other cosmetic changes that do not interfere with the designed function of the unit.

This warranty does not include failure of the product due to alteration, modifications, misuse, abuse, improper maintenance, improper installation or repair.

This warranty is not transferable and is only applicable to the original purchaser named on a completed and registered warranty record at JK-Products, Inc. headquarters.

JK Products, Inc. obligations do not include transportation charged on the shipping of new components to the customer or on the return of defective components of the manufacturer. All freight charges must be paid by the customer.

There are no warranties which extend beyond the description on the face hereof. There is no express or implied warranty of fitness or merchantability.

The remedies provided in this limited warranty are the exclusive remedies provided to the purchaser by Global Service USA and are provided in substitution of all other remedies. Consequential and incidental damages are excluded.



NOTE:

Any misuse, abuse, modification, alteration, improper maintenance, improper installation, etc. of this product or its intended use will subject said party to direct responsibility and liabilities for any and all consequences resulting from such action. Furthermore, such actions will void all warranties and may relieve Global Service USA of any liabilities relating to said product.

18 Index



A

Acrylic panel.....	73
Load	30
Air conditioning	
Filter	90
Air intake openings	95
Air routing	16
Application.....	56
Application tips	61
AQUA / AROMA SYSTEM maintenance.....	93
AQUA SYSTEM canister	94
AROMA canister	93

B

Bed base	76
Bluetooth®	58

C

Canopy panel.....	74
Changing the chip cards.....	78
Cleaning	62, 66
Air conditioning.....	90
Filter	88
Filter mats.....	89
Lamps, filter panels.....	67
Cleaning or replacing lamps	
Base.....	82
Cleaning or replacing the lamps	
Shoulder tanner.....	86
Cleaning/disinfection	
Cleaning intervals	68
Commissioning	18
Condensation tank.....	90
Control Center	48
Customer Service	2

D

Decommissioning	96
Definitions.....	7
Device description	13
Dimensions.....	32
Disinfection.....	62, 66
Disposal.....	96



E

Electrical specifications	30
EMC information	102
Export	12

F

Facial tanner	
Cleaning/replacing	84
Maintenance intervals	69
Filter	
Cleaning	88, 89
Filter panels	81, 83
Cleaning/replacing	84
Cleaning/replacing	86
Maintenance intervals	69
Foreseeable misuse	41
Functions	51

G

Guidelines	12
------------------	----

H

High pressure lamps, fire hazard	81
--	----

I

Installation location	16
Intended use	35

L

Lamp configuration	33
Lamps	
Control	33
Power rating	33
Side section	83
Lamps cleaning/replacing	
Facial tanner	84
Lamps, replacing	77
Low pressure lamps	79

N

Noise level	30
-------------------	----



O

Obligations of the studio operator.....	41
Operating altitude	17
Operating display of the sensor	57
Order numbers.....	33, 34

P

Packaging.....	96
Plastic surfaces (care).....	67
Power rating.....	30
Prohibited usage.....	36

R

Recycling	96
-----------------	----

S

Safety measures during maintenance work	63
Scope of delivery	12
Service kit.....	72
Setup location	17
Side section panel	75
Skin condition.....	35
Skin types.....	46
Spare parts.....	69
Order numbers.....	34
Staff training	41
Starter	33, 80
Storage	96

T

Tanning times.....	47
Technical data.....	30
Transport.....	15
Troubleshooting	99

V

Ventilation.....	16
------------------	----

